


Trident Integral HPLC Guard Column System

The Trident system's foundation consists of the analytical column configured with our exclusive Trident end fitting and XF end fitting. You can order any Restek HPLC column with this configuration simply by adding a "-700" suffix to the catalog number for the column. The external cap frit in the XF end fitting can be replaced by removing the XF end fitting, pulling off the old frit, and pushing a new frit into place. When replacing the XF end fitting, tighten less than 1/8 turn past finger-tight—there is no need to overtighten this connection.

For maximum protection against particulate matter and sample contaminants, the system can be configured with both an integral guard cartridge and a replaceable external frit. To assemble this configuration, remove the XF end fitting, insert the appropriate 1 cm or 2 cm guard cartridge into the Trident end fitting, attach an XG-XF fitting (cat.# 25026 or 25062) to the Trident end fitting, and connect the XF end fitting to the XG-XF fitting. Do not overtighten any of these connections.



Assembled column with Trident integral guard system. To order, add "-700" to the catalog number of the column, and order guard cartridges and the XG-XF fitting separately.

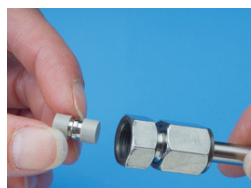

Remove the XF end fitting and install the guard cartridge in the end of the column.



Add the XG-XF fitting.



Reinstall the XF end fitting with cap frit.



The cap frit can be easily replaced if it becomes contaminated/plugged.

Replacement parts for the Trident guard column system

Part	cat.#
XF End Fitting	25024
XG-XF Fitting for 10 mm Guard Cartridge	25026
XG-XF Fitting for 20 mm Guard Cartridge	25062
Cap Frits, 4 mm/2.0 µm porosity	25022
Cap Frits, 4 mm/0.5 µm porosity	25023
Cap Frits, 2 mm/2.0 µm porosity	25057

Restek guard column cartridges also can be used with Trident Direct holders:

Part	cat.#
10 mm Guard Cartridge Holder with Filter	25084
20 mm Guard Cartridge Holder with Filter	25086

Contact Technical Service at 1-800-356-1688, 1-814-353-1300, ext. 4, or support@restek.com (or contact your Restek representative) if you have any questions about this product or any other Restek product.



Restek U.S. • 110 Benner Circle • Bellefonte, PA 16823
 phone: 1-814-353-1300 • 1-800-356-1688 • fax: 1-814-353-1309 • www.restek.com
Restek France • phone: +33 (0)1 60 78 32 10 • fax: +33 (0)1 60 78 70 90 • www.restek.fr
Restek GmbH • phone: +49 (0)6172 2797 0 • fax: +49 (0)6172 2797 77 • www.restekgmbh.de
Restek Ireland • phone: +44 (0)2890 814576 • fax: +44 (0)2890 814576 • e-mail: restekireland@aol.com
Restek Japan • phone: +81 (3)6459 0025 • fax: +81 (3)6459 0025 • e-mail: restekjapan@restek.com
Thames Restek U.K. LTD • phone: +44 (0)1494 563377 • fax: +44 (0)1494 564990 • www.thamesrestek.co.uk