


Liquid Syringes

Series CG-130

Fixed needle: .019" x .005" x 2", bevel, open end
 Removable needle: .019" x .005" x 2.25", bevel, open end

- ▶ **Reinforced plunger** – eliminates plunger bending
- ▶ **Plunger guide in barrel** – smaller front section of plunger is confined within guide section
- ▶ **Replaceable needles** – low dead volume connection for a leak-tight fit
- ▶ **Leak-tight to 250 psi** – liquids and gases

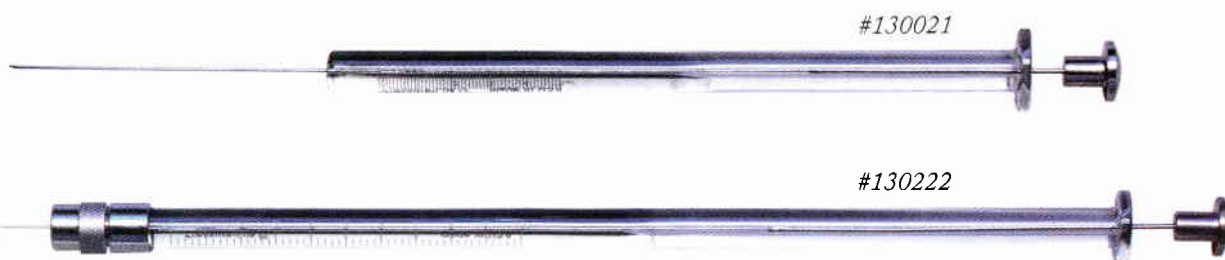

A heavy, reinforced plunger within an enlarged plunger guide section of the CG-130 barrel virtually eliminates plunger bending. The smaller front section of the plunger is confined within the guide section even when the plunger is retracted, while the sturdy upper section absorbs the force of a rapid downstroke.

	5 µl	10 µl	Replacement needles
Fixed needle:	130021	130022	Size: .019" x .005" x 2.25"
			Type: Bevel, open end
Removable needle:	130221	130222	Pkg. of 3: 123050

With the Series CG-130
 you'll never bend another plunger