

A




B




C



D

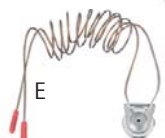


E

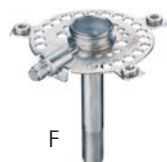

F



G

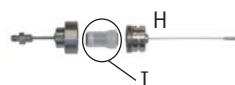

H

I

Direct Replacement and Siltek® Treated Split/Splitless Injection Ports for Agilent GCs

Would you like better performance from your injector? Our Siltek® treated split/splitless injector is a direct replacement for Agilent 5890 and 6890/6850 GCs, but Siltek® treatment passivates the metal surface to ensure an inert pathway for the sample, delivering increased performance. The injector is manufactured from high-quality stainless steel and meets or exceeds Agilent original equipment performance.

Direct Replacement and Siltek® Treated Split/Splitless Injection Port for Agilent 5890 GCs

Description	Similar to Agilent part #	qty.	cat.#
A) Replacement Weldment*	19251-60575	ea.	20265
Siltek® Weldment**	19251-60575	ea.	20267
B) Replacement Shell Weldment	19251-80570	ea.	20266
Siltek® Shell Weldment**	19251-80570	ea.	20268
C) O-rings for Agilent Trap Fittings	5180-4181	25-pk.	22064

Direct Replacement and Siltek® Treated Split/Splitless Injection Port for Agilent 6890/6850 GCs

Description	Similar to Agilent part #	qty.	cat.#
D) Replacement Weldment with EPC	G1544-60575	ea.	22674
Siltek® Weldment with EPC**	G1544-60575	ea.	22672
E) Replacement Weldment*	19251-60575	ea.	20265
Siltek® Weldment**	19251-60575	ea.	20267
F) Replacement Shell Weldment	G1544-80570	ea.	22673
Siltek® Shell Weldment**	G1544-80570	ea.	22671
G) Optional Split/Splitless Weldment (for use with large canister type filter)	G1544-60585	ea.	22686
Siltek® Optional Split/Splitless Weldment** (for use with large canister type filter)	G1544-60585	ea.	22670

*For use with manual flow or EPC on Agilent 5890 GCs. For use with manual flow only on Agilent 6890/6850 GCs.

**Treated with Restek's exclusive Siltek® treatment for increased inertness.

Replacement Chemical Traps for Agilent GCs

- Easy to install.
- Attach to same fittings as original equipment.
- Built-in frits retain fine particles; adsorbents remove both moisture and hydrocarbons.

Description	Similar to Agilent part #	qty.	cat.#
H) Optional Split Vent Trap Assembly for Agilent 6890/6850 GCs	G1544-60610	kit	23031
I) Replacement Traps (2) and O-rings (4)	G1544-80530	kit	23032

also available

Looking for autosampler syringes?

Restek carries a full line of autosampler syringes for Agilent instruments from both Hamilton and SGE. See [pages 278 to 280](#) or visit www.restek.com and search on syringes.

