


Organochlorine Pesticides
US EPA Method 8081
Rtx®-35

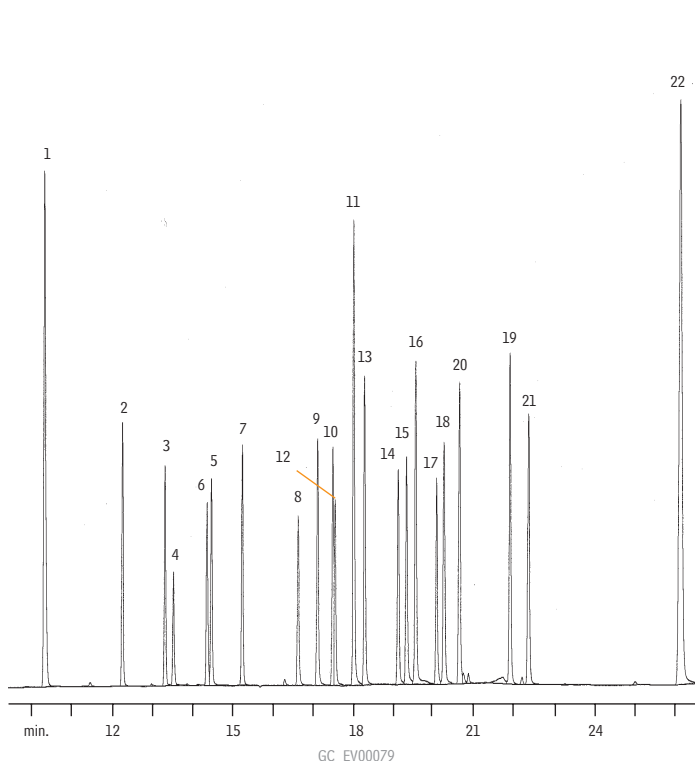

1. 2,4,5,6-tetrachloro-*m*-xylene (SS)
2. α -BHC
3. γ -BHC (lindane)
4. β -BHC
5. δ -BHC
6. heptachlor
7. aldrin
8. heptachlor epoxide
9. γ -chlordane
10. α -chlordane
11. 4,4'-DDE
12. endosulfan I
13. dieldrin
14. endrin
15. 4,4'-DDD
16. endosulfan II
17. 4,4'-DDT
18. endrin aldehyde
19. methoxychlor
20. endosulfan sulfate
21. endrin ketone
22. decachlorobiphenyl (SS)

GC_EV00079

Column: Rtx®-35, 30m, 0.32mm ID, 0.50 μ m (cat.# 10439)
Oven temp.: 120°C (hold 1 min.) to 285°C
@ 8.5°C/min. (hold 6 min.)
Inj.: direct, Uniliner® inlet liner
(cat.# 20335), 200°C
Det.: ECD, 300°C with anode purge
Flow rate: helium @ 2.1mL/min. set @ 120°C



Catch the Buzz!

To automatically receive free technical literature electronically, sign up for Restek's popular e-newsletter, *The Buzz*, at www.restek.com/buzz