



**Aroclor 1248 PCBs
Rtx[®]-CLPesticides & Rtx[®]-CLPesticides2**

Rtx[®]-CLPesticides

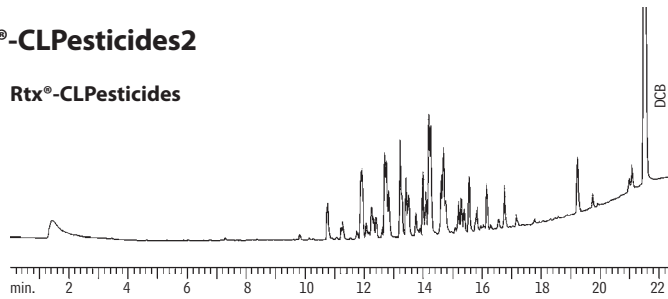

Rtx[®]-CLPesticides2

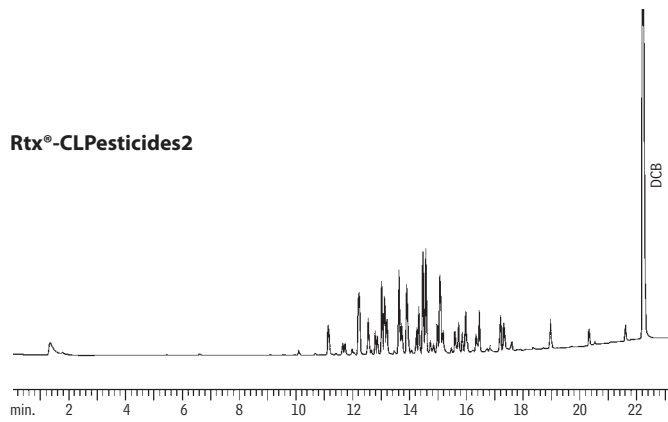


GC_EV00453

Columns: Rtx[®]-CLPesticides, 30m, 0.32mm ID, 0.50µm (cat.# 11139)
 Rtx[®]-CLPesticides2, 0.32mm ID, 0.25µm (cat.# 11324)
 Inj.: 0.5µL direct injection, 160pg on-column concentration
 Oven temp.: 120°C (hold 1 min.) to 300°C @ 9°C/min. (hold 10 min.)
 Inj. temp.: 210°C
 Det.: ECD with anode purge, 310°C
 Carrier gas: helium

**Aroclor 1254 PCBs
Rtx[®]-CLPesticides & Rtx[®]-CLPesticides2**

Rtx[®]-CLPesticides

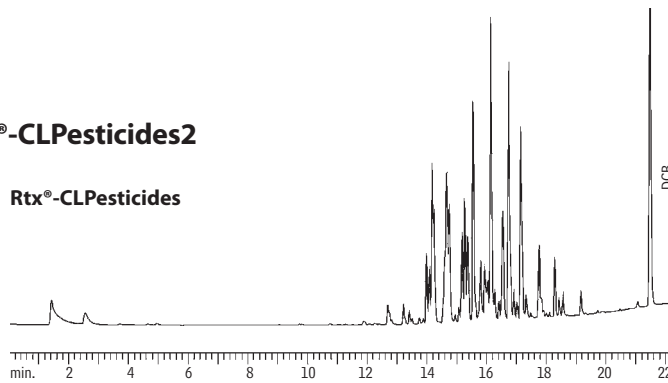

Rtx[®]-CLPesticides2

