


Liner-to-Liner Reproducibility: A Measure of Consistent Quality

It's not good enough to have one quality liner. You have to be confident that every liner will give the same level of performance. We test Sky™ liners extensively to ensure that each one is exceptionally inert and will provide optimal results. Using endrin breakdown as a measure of reproducibility, the data in Figure 4, based on multiple lots, illustrate that Sky™ liners are more consistently inert than competitor products.

New Sky™ liners provide exceptional inertness across a wide range of active analytes. The consistent, comprehensive deactivation process results in the accuracy and precision you need for reliable trace level analyses. Simplify liner selection with Sky™ liners from Restek—**choose blue, the best choice for dependable results.**

Figure 4 Sky™ inlet liners from Restek consistently show less endrin breakdown than comparable liners from other sources.

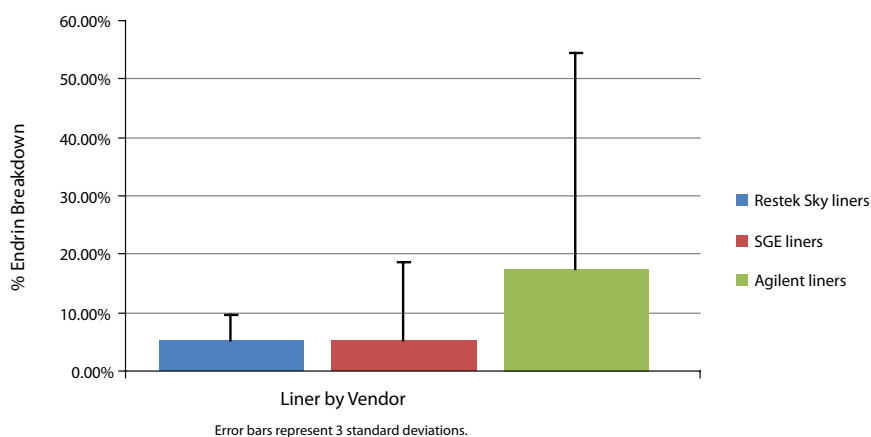

New Sky™ inlet liners are easily recognizable for your reordering convenience. All Sky™ liners come in specially marked boxes and are packaged in ultra-clean blister packs.