

**Column Hanger**  
for Agilent 7890, 6890, 5890,  
and 5880A GCs




Description	qty.	cat. #	price
Column Hanger for Agilent 7890, 6890, 5890, and 5880A GCs	ea.	22128	

**Sapphire Scribe**

- Cuts fused silica tubing.
- Produces a clean, square cut.



20182

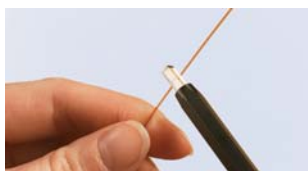

One quick stroke...



...and just a tap leaves a clean, square end.

Description	qty.	cat.#	price
Sapphire Scribe	ea.	20182	

**Scoring Wafer with Handle**

- Ceramic wafer is serrated on one side and straight-edged on the other to cut both fused silica and metal tubing cleanly.
- Unique, ergonomic handle is made of soft, comfortable rubber.



23015

Description	qty.	cat.#	price
Scoring Wafer with Handle	2-pk.	23015	



23027



23026

**Shortix® Capillary GC Column Cutter**

- Consistently make precise, clean, square cuts with a diamond blade.
- Built-in magnifier to verify square cut.
- Use with 0.25 mm ID to 0.53 mm ID tubing (0.78 mm OD maximum).


Description	qty.	cat.#	price
Shortix Capillary GC Column Cutter	ea.	23026	
Maintenance Kit for Shortix Capillary GC Column Cutter (Includes: diamond cutting wheel, O-rings, and a tool to open the column cutter)	kit	23027	

**Ceramic Scoring Wafer**

- Four straight scoring edges for cutting fused silica tubing and four serrated edges for cutting MXT® metal capillary columns.
- Sure-grip handle included.



20116

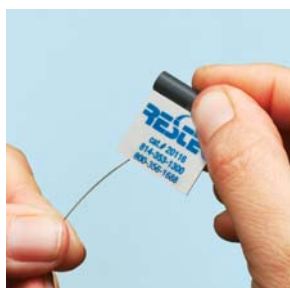


Exert just enough pressure to put a slight arc in the tubing. The tubing should fall off or break with a slight tap of the wafer.



Check the cut against the white of the scoring wafer. Look for a clean, square cut.

Description	qty.	cat.#	price
Ceramic Scoring Wafers	5-pk.	20116	

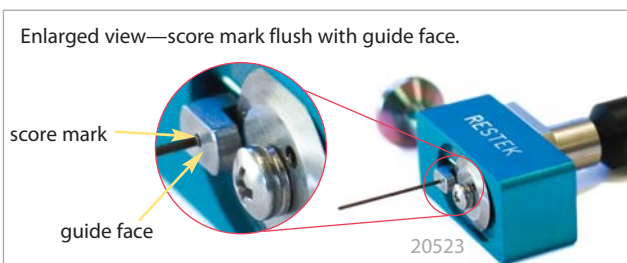

**Restek Tubing Scorer**

- Makes perfect cuts every time.
- Easy to use.
- Leaves column entrance perfectly round.

Metal MXT® columns are easy to cut. Scoring wafers can be used, but may leave the column end irregularly shaped. The Restek tubing scorer is designed to make a perfect cut every time, leaving the column entrance perfectly round.

Description	qty.	cat.#	price
Restek Tubing Scorer for MXT Columns (0.25-0.53mm ID & 0.5-0.8mm OD)	ea.	20523	
Replacement Scoring Wheel	ea.	20522	

