


restek
innovation!



Fit both Restek cage designs.

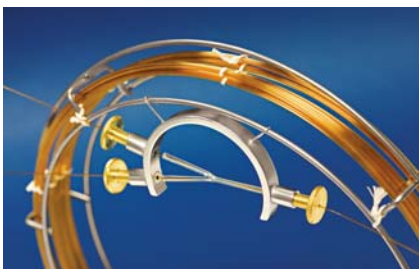

Patent pending.



SeCure® “Y” Connector Kits

- Connect two analytical columns to a transfer line or guard column.
- Use standard “Y” Press-Tight® connectors and 1/16" graphite ferrules.
- Reliable seal integrity, will not unexpectedly disconnect during temperature-programmed analyses.
- Open design allows visual confirmation of the seal for added confidence in the connection.
- Fit both Restek cage designs.

Combine the simplicity of a “Y” Press-Tight® connector with the strength of a metal union. The ferrules and knurled nuts hold the fused silica tubing in place, which prevents the tubing from unexpectedly disconnecting, even at temperatures as high as 400 °C.

The SeCure® “Y” Connector’s open design allows visual confirmation of the seal.

Dual-column confirmational analysis with a single injection—one of the SeCure® “Y” connector’s many uses.

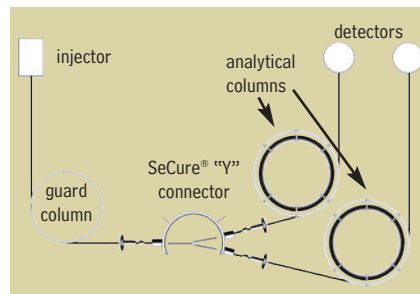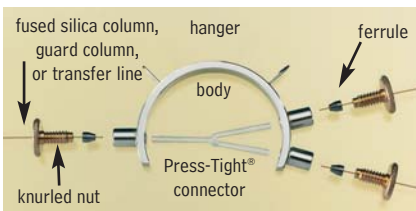

Secure, reliable
column-to-column
connections!


Kits include: SeCure® “Y” connector body, 3 knurled nuts, “Y” Universal Press-Tight® union, 3 ferrules.

Description	Ferrules Fit Column ID	qty.	cat.#	price
SeCure “Y” Connector Kit	0.18/0.25/0.28mm	kit	20276	
SeCure “Y” Connector Kit	0.32mm	kit	20277	
SeCure “Y” Connector Kit	0.45/0.53mm	kit	20278	
Knurled nut		3-pk.	20279	

Universal “Y” Press-Tight® Connectors

An alternative method of performing dual-column confirmational analyses!

- Split sample flow onto two columns.
- Split a single column flow to two detectors—perform confirmation analysis with a single injection.
- Deactivated Press-Tight® connectors assure better recovery of polar and nonpolar compounds.
- Siltek® treated connectors are ideal for organochlorine pesticides analysis.
- Fit column ODs from 0.33–0.74 mm (Restek 0.1 mm–0.53 mm ID).

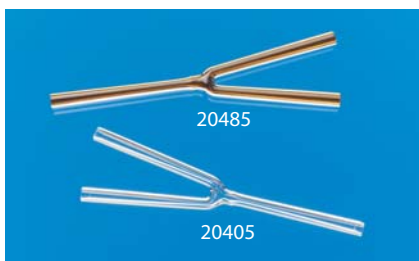

Description	ea./price	3-pk./price
Universal “Y” Press-Tight Connector	20405 \$78	20406
Universal “Y” Press-Tight Connector, Deactivated	20405-261 \$79	20406-261
Universal “Y” Press-Tight Connector, Siltek Treated	20485 \$80.50	20486

Graphite Ferrules

for SeCure® “Y” Connectors

- Preconditioned to minimize out-gassing.
- High-purity, high-density graphite.
- Stable to 450 °C.
- No binders that can off-gas or adsorb analytes.
- Smooth surface and clean edges.



Ferrule ID	Fits Column ID	Graphite 10-pk./price	Graphite 50-pk./price
0.4mm	0.10/0.15/0.18/0.25/0.28mm	20200 \$32	20227
0.5mm	0.32mm	20201 \$32	20228
0.8mm	0.45/0.53mm	20202 \$32	20224