

TUBING

Pre-cut stainless steel tubing

UHPLC

SPECS

- Pressure rating varies with ID; minimum is 365 bar (5300 psi) (see page 93)
- Tolerances
 - OD ±0.05mm (.002")
 - ID ±0.025mm (.001")

- For zero dead volume connections
- Steam-cleaned and ready to use
- Soft annealed OD for easy ferrule connections
- Color bands for quick identification of tubing ID

VICI Jour offers pre-cut 1/16" stainless steel tubing in lengths most commonly required in UHPLC and HPLC systems. The tubes are cut to square, burr-free ends for zero dead volume connections. Colored bands allow quick ID identification.

VICI pre-cut tubing is specially cleaned with steam from micro-filtered deionized water, which removes both organic and inorganic contaminants. Longer tubes are cleaned with liquid hexane. Softener-free polyethylene caps are used to prevent any contamination during shipment and storage.

Length:	50 mm		100 mm		200 mm		300 mm	
OD Color	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	
1/16"	Red	0.13	JR-T-97005	JR-T-97010	JR-T-97015	JR-T-97020	JR-T-97025	
	Yellow	0.18	JR-T-97025	JR-T-97030	JR-T-97035	JR-T-97040	JR-T-97045	
	Blue	0.25	JR-T-97045	JR-T-97055	JR-T-97060	JR-T-97065	JR-T-97070	
	Orange	0.50	JR-T-97075	JR-T-97080	JR-T-97085	JR-T-97090	JR-T-97095	
	Green	0.75	JR-T-97095	JR-T-97100	JR-T-97105	JR-T-97110	JR-T-97115	



Pre-cut flexible stainless steel tubing

UHPLC

SPECS

- Pressure rating varies with ID; minimum is 365 bar (5300 psi) (see page 93)
- Tolerances
 - OD ±0.05mm (.002")
 - ID ±0.025mm (.001")

- Seamless 1/32" tubing with 1/16" laser-welded ends
- For zero dead volume connections
- Steam-cleaned and ready to use
- Soft annealed OD for easy ferrule connections
- Color bands for quick identification


Our flexible pre-cut stainless steel tubing is made of 1/32" premium grade 316 stainless steel tubing with ends laser-welded into 1/16" tubing of the same premium grade. Available in lengths most commonly required in UHPLC and HPLC systems.



The tubes are cut to square, burr-free ends for zero dead volume connections. Colored bands allow quick ID identification.

VICI pre-cut tubing is specially cleaned with steam from micro-filtered deionized water, which removes both organic and inorganic contaminants. Longer tubes that cannot be steam-cleaned are cleaned with liquid hexane. Softener-free polyethylene caps are used to prevent any contamination during shipment and storage.

Length:	105 mm		150 mm		200 mm		280 mm		400 mm	
OD Color	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	
1/16"	Red	0.13	JR-T-96005	JR-T-96010	JR-T-96015	JR-T-96020	JR-T-96025	JR-T-96030	JR-T-96035	
	Yellow	0.18	JR-T-96105	JR-T-96110	JR-T-96115	JR-T-96120	JR-T-96125	JR-T-96130	JR-T-96135	
	Blue	0.25	JR-T-96205	JR-T-96210	JR-T-96215	JR-T-96220	JR-T-96225	JR-T-96230	JR-T-96235	
	Orange	0.50	JR-T-96305	JR-T-96310	JR-T-96315	JR-T-96320	JR-T-96325	JR-T-96330	JR-T-96335	
		600 mm		700 mm		800 mm		900 mm		
		Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	Product No.	
	Red	0.13	JR-T-96030	JR-T-96035	JR-T-96040	—	—	—	—	
	Yellow	0.18	JR-T-96130	JR-T-96135	JR-T-96140	JR-T-96145	—	—	—	
	Blue	0.25	JR-T-96230	JR-T-96235	JR-T-96240	—	—	—	—	
	Orange	0.50	JR-T-96330	JR-T-96335	JR-T-96340	—	—	—	—	

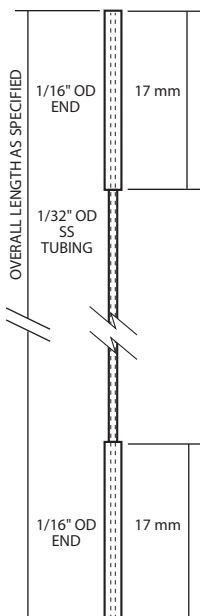

TECH TIP

Contact your local distributor or VICI for custom lengths.