



## Tools

### Custom socket wrench

This 1/4" socket wrench with a slot to slip over 1/16" tubing works great for all types of 1/4" hex nuts (such as Valco 1/16" ZDV fitting nuts). It's especially useful when nuts are difficult to access with an open end wrench.

Prod No      Price  
SWH4



### Ferrule removal kit

When polymeric ferrules get stuck in a fitting detail, these little ferrule spears will save you from becoming so irritated that you tear up your entire lab in frustration. Each kit includes two sizes of tapered stabbers for retrieving capillary size ferrules.

Prod No      Price  
FRK1



### Hex key set

The hex key set has a wrench to fit any socket head screw on any VICI valve or actuator. Includes the following sizes: .050", 1/16", 5/64", 3/32", 7/64", 1/8", 9/64", and 5/32".

Prod No      Price  
HKS



#### TECH TIP

If a fused silica tube breaks off in a through-type union, remove the nuts and the tube opposite the broken one. Clear the fitting by passing a drill or wire of the appropriate diameter into the unbroken side and through the center of the fitting.

Our ferrule removal kit can be used to remove ferrules from tee and cross fittings.




2009 #60



**Open end wrenches**

Size	For use with	Prod No	Price
3/16" x 1/4"	1/32" and 1/16" nuts	OEW	
3/8" x 7/16"	1/8" nuts	OEW-2	
1/2" x 9/16"	1/4" nuts	OEW-3	

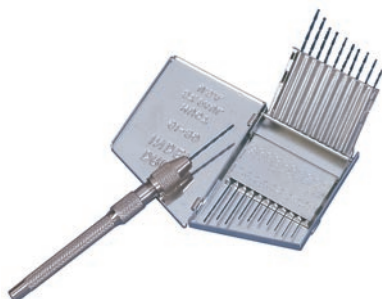

**Pin vise and drill index**


The drill index has drills sized from 0.0135" to 0.039" (0.34 to 1 mm). These are useful tools when a fused silica tube breaks in a union (see Tech Tip on the facing page) and for enlarging the inner diameter of fused silica adapters.

Prod No	Price
PV	

**Template**

This tool is just what you need when you're working out plumbing and valve switching schematics. It features templates for two position valves with 4, 6, 8, and 10 ports with indications of both positions, as well as various flow symbols. For added convenience, the sides are edged with metric and inch rulers.

Prod No	Price
TEMPLATE1	

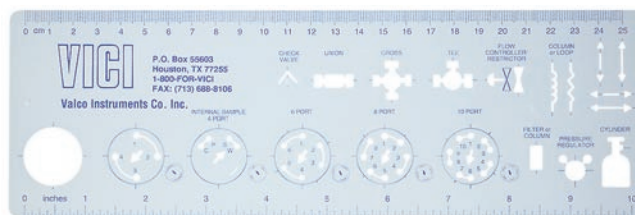

**MORE INFORMATION**

- Tools for valves
- Pencil magnet . . . . p 210
- Valve spanner handle. . . . . 211
- Tightening tools for PEEK fittings. . . . . 67
- Tubing accessories . . . . 90

CHROMalytic TECHNOlogy Pty Ltd AUSTRALIAN Distributors e-mail: sales@chromtech.net.au Tel: 03 9762 2034