




Flanging Tool, Starter Kit

Cheminert flanging tools

The flanging tool makes the flange which retains the standard 1/4-28 tube end fitting and washer on PTFE or FEP tubing. With this tool, lengths of tubing may be easily assembled to any required dimension. The time required is approximately 5 to 10 seconds per flange.

Flanging tools are available for 110 VAC or 230 VAC, and come complete with tips for 0.75 mm, 1.0 mm, and 2.00 mm ID tubing, a tubing holder for gripping the tubing during the flanging operation, a razor blade for tube cutting, and instructions.

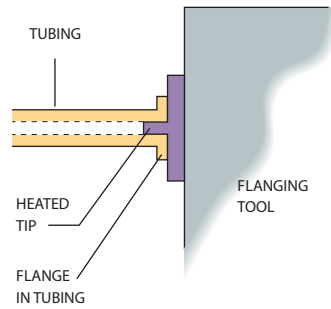

Flange being made on tubing

	<i>Prod No</i>	<i>Price</i>
Flanging tools 110 VAC	CFT-110	
230 VAC	CFT-220	
Flanging tool accessories		
Flanging tips		
for tubing ID ≤ 0.25 mm	CFT-TXC	
for tubing ID ≤ 0.75 mm	CFT-TC	
for tubing ID ≤ 1.00 mm	CFT-TM	
for tubing ID ≤ 1.50 mm	CFT-TL	
for tubing ID ≤ 2.00 mm	CFT-TXL	
Razor blades (<i>pkg /10</i>)	CFT-R	
Tubing holder	CFT-H	



Cheminert starter kits

Starter kits come in either 1/16" or 1/8" versions, with flanging tools for 110 VAC or 230 VAC.

	110 VAC		230 VAC	
	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>
Starter kits				
1/16" tubing	CFT1K-110		CFT1K-220	
1/8" tubing	CFT2K-110		CFT2K-220	

The starter kit includes:

- 1 flanging tool with 3 flanging tips
- 1 tubing holder
- 20 standard tube end fittings
- 20 stainless steel washers
- 10 couplings
- 20 feet of PTFE tubing
(1/16" OD x 0.030" ID
or
1/8" OD x .060" ID)
- 1 male luer adapter
- 1 female luer adapter
- 1 plug
- 1 tee
- 1 glass connector



MORE INFORMATION

Standard tube end fittings	page 69
Stainless steel washers	69
Couplings	72
Male luer adapter	76
Female luer adapter	76
Plug	71
Tee	74
Glass connector	75