


Polymeric Tubing

Polymeric tubing is square cut and ready to use. Each package of polymeric tubing contains one piece of the specified length.

See also PEEK tubing, pages 88-89.

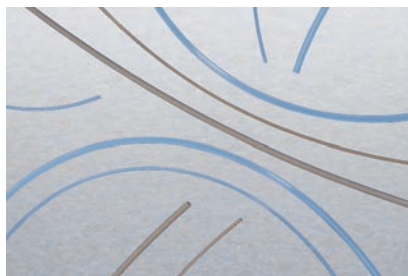

1/16" OD polymeric tubing

	.006" ID		.010" ID		.015" ID		.020" ID		.030" ID	
	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>
PTFE										
10 feet	TTF106-10FT		TTF110-10FT	\$15	TTF115-10FT		TTF120-10FT		TTF130-10FT	
25 feet	TTF106-25FT		TTF110-25FT	35	TTF115-25FT		TTF120-25FT		TTF130-25FT	
100 feet	TTF106-100FT		TTF110-100FT	130	TTF115-100FT		TTF120-100FT		TTF130-100FT	
FEP										
10 feet	-		TFEP110-10FT	16	-		TFEP120-10FT		TFEP130-10FT	
25 feet	-		TFEP110-25FT	39	-		TFEP120-25FT		TFEP130-25FT	
100 feet	-		TFEP110-100FT	156	-		TFEP120-100FT		TFEP130-100FT	
ETFE										
10 feet	-		TTZ110-10FT	25	-		TTZ120-10FT		TTZ130-10FT	
25 feet	-		TTZ110-25FT	57	-		TTZ120-25FT		TTZ130-25FT	
100 feet	-		TTZ110-100FT	210	-		TTZ120-100FT		TTZ130-100FT	

1/8" OD polymeric tubing

	.030" ID		.060" ID		.085" ID	
	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>	<i>Prod No</i>	<i>Price</i>
PTFE						
10 feet	TTF230-10FT		TTF260-10FT		TTF285-10FT	
25 feet	TTF230-25FT		TTF260-25FT		TTF285-25FT	
100 feet	TTF230-100FT		TTF260-100FT		TTF285-100FT	
FEP						
10 feet	-		TFEP260-10FT		-	
25 feet	-		TFEP260-25FT		-	
100 feet	-		TFEP260-100FT		-	
ETFE						
10 feet	-		TTZ260-10FT		-	
25 feet	-		TTZ260-25FT		-	
100 feet	-		TTZ260-100FT		-	

- Not available


MORE INFORMATION

PEEK tubing
 Natural page 88
 Color-coded 89
 Striped 89

CLEAN-CUT POLYMER TUBING CUTTER

For leak-free tubing connections in an LC system, right angles and clean cuts are essential. The Clean-Cut makes burr-free perpendicular cuts on polymeric tubing without distorting the outside diameter or closing the inside diameter. The handy pocket-sized tool features a unique safety locking mechanism to secure the blade when not in use.

	<i>Prod No</i>
Clean-Cut tubing cutter	JR-797
Replacement blade	JR-798

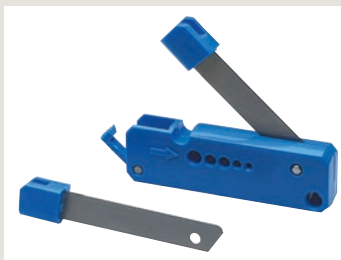

TUBING CLIP – THE LC TUBING ORGANIZER

The tubing clip holds 1/16" and 1/8" polymer tubing precisely where you want them in your beakers, flasks, bottles, etc. up to 4 mm wall thickness. The stainless steel spring ensures a long lifetime.

Package of 5:	<i>Prod No</i>
Tubing clip	JR-9001-5

TUBING POLYMERS

PTFE Inert; very soft, easily cold flows. Produced as Teflon®

FEP Chemically resistant like PTFE, but lower creep and higher friction. More transparent than PTFE.

ETFE Resistant to most chemical attack; some chlorinated solvents will cause tubing to swell. Produced as Tefzel®

10 ft	=	3.05 m
25 ft	=	7.62 m
100 ft	=	30.48 m