


Selectors – Low Pressure

Trapping/flow-through flowpath – STF configuration

The STF selector is a variation of the ST flowpath, with the single difference that the non-selected streams are returned to their own vents or sources rather than being dead-ended or trapped as they are in the standard ST configuration. For an application suggestion, see page 138.

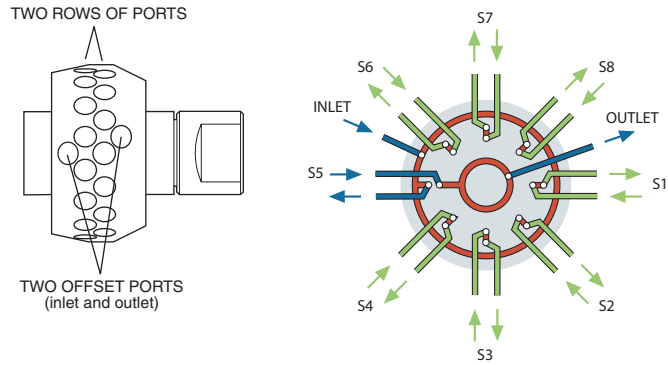

1/16" fittings, 0.75 mm ports (.030")

MW Type

Low pressure	Includes 2" standoff. Ask about closemount assembly if valve will not be heated.	Standard electric actuators: 110 VAC for USA 110/230 VAC to 24 VDC power supply for international	SPECS 200 psi gas 200°C max Nitronic 60 body Valcon E rotor
STF Trap/ flow-throw		Microelectric actuators: 24 VDC (includes a 110/230 VAC to 24 VDC power supply)	
1/16" 0.75 mm			

OPTIONS

- 4 and 8 positions available
- 3", 4", and 6" standoffs
- Materials: Hastelloy C, Inconel 600, Monel 400, Nickel 200, Nitronic 50, Titanium (see pages 254-255)

	6 Position		10 Position		12 Position		16 Position	
	Prod No	Price	Prod No	Price	Prod No	Price	Prod No	Price
Manual (not recommended)	2CSTF6MWE		2CSTF10MWE		2CSTF12MWE		2CSTF16MWE	
With air actuator	A2CSTF6MWE		A2CSTF10MWE		A2CSTF12MWE		A2CSTF16MWE	
With standard elec actuator	E2CSTF6MWE		E2CSTF10MWE		E2CSTF12MWE		E2CSTF16MWE	
With microelectric actuator	EMT2CSTF6MWE		EMT2CSTF10MWE		EMT2CSTF12MWE		EMT2CSTF16MWE	
Replacement valve	DCSTF6MWE		DCSTF10MWE		DCSTF12MWE		DCSTF16MWE	
Replacement rotor	SSACSTF6MWE		SSACSTF10MWE		SSACSTF12MWE		SSACSTF16MWE	

MORE INFORMATION

- Application. page 138
- Actuators
 - Air 194
 - Microelectric .. 190-191
 - Standard elec. 193
- Materials
 - Metals. 254-255
 - Polymers 256
 - Valve rotors. 257
- Mounting hardware
 - Closemount 208
 - Standoff. 205

CHROMalytic TECHnology Pty Ltd AUSTRALIAN Distributors e-mail: sales@chromtech.net.au Tel: 03 9762 2034