



Knobs, Handles, and Right Angle Drives

Knobs and handles

for use with a standoff

If you already have a spare standoff assembly (see facing page) but lack the knob or retainer, or have an actuated valve on a standoff which you'd like to convert to manual use, here's what you'll need. Includes knob or handle, retainer, and adapter.

Knob for a W type valve	Prod No WMMASO
T-handle for a UW type valve	SOUTH

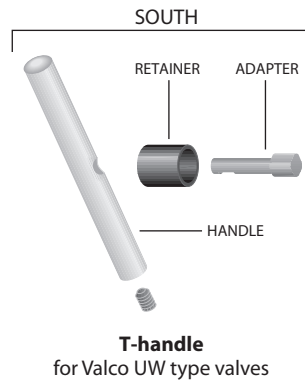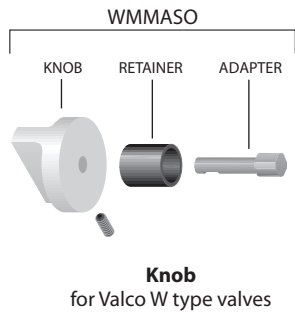

RAD Right angle drive

for two position actuators

Some installations don't allow the valve and actuator to be installed in a typical in-line configuration. The RAD is a 90° gearbox which permits the actuator or handle to be installed at a right angle to the valve. The RAD fits all VICI two position electric and air actuators. Not for use with 1/4" valves.

RAD with standoff includes a 2" standoff; 3", 4", and 6" standoffs are also available.

With closemount hardware	With 2" standoff assembly
Prod No RAD	Prod No 2RAD

