




VALUE PLASTICS, INC.
SPECIALISTS IN FLUID CONNECTIONS™

THE SF
SERIES
SANITARY FITTINGS

SF Series sanitary fittings from Value Plastics are designed for use in single-use container fluid transport systems, manifolds and wherever sensitive fluids must be transferred. They are offered in two flange sizes, 1-1/2" ("Maxi"), and the popular 3/4" ("Mini") size commonly used with small tubing sizes, both work with flanges made to the ISO 2852 and ASME-BFE specifications*. SF Series sanitary fittings incorporate 600 Series barbs for easy use with silicone, C-Flex® or Tygon® tubing.

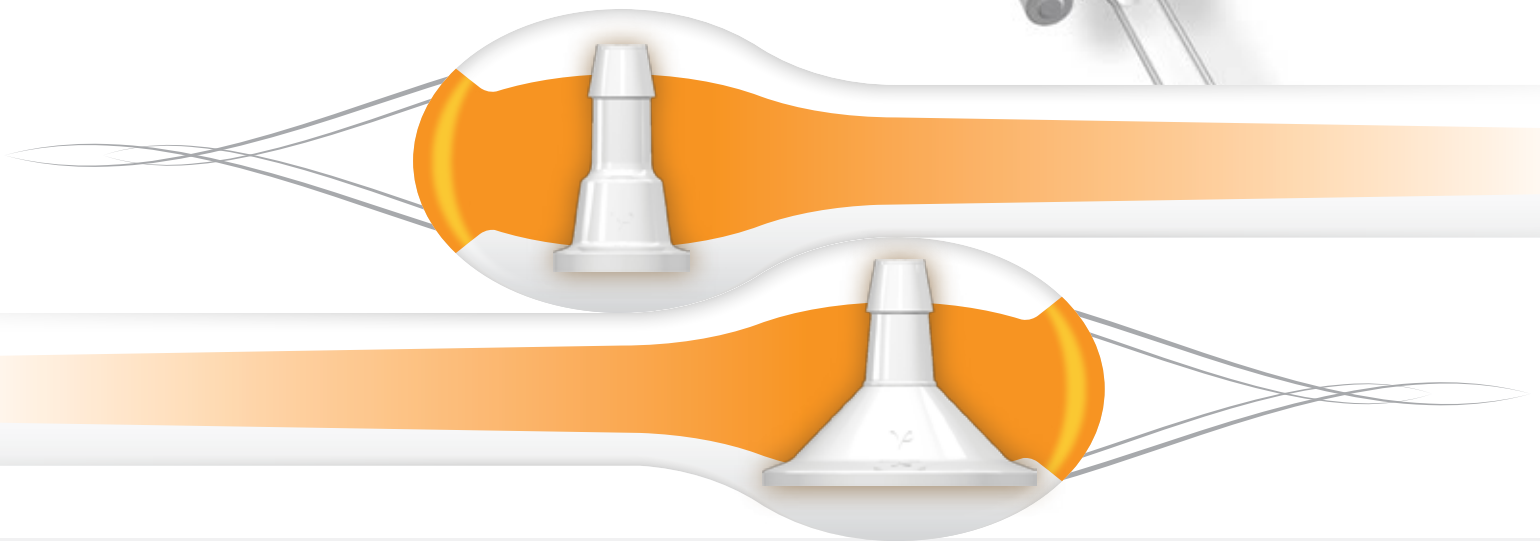

FEATURES

BENEFITS

SPECIFICATIONS

Animal-free materials	No BSE documentation required	Materials: Animal Free Polypropylene, Natural PVDF Tubing Sizes: 1/4", 3/8", 1/2", 3/4", 1"
Radiused flow surfaces	Smooth uninterrupted flow	Leak Resistance: Tested at 30 psi air under water for 90 seconds with no leaks Tubing Removal Force: 1/4" PP with C-Flex shore tubing: 21 lbs. 1/2" PP with C-Flex shore A tubing: 89 lbs. 1" PP with C-Flex shore A tubing: 121 lbs.
"Mini" & "Maxi" sizes	Easy to retrofit current operations Works with commercially available clamps	Interior Surface Finish: 15 Ra Microinch
No Particulate build-up	Interior surface controlled to 15 Ra microinch	

For more information and samples, please visit www.valueplastics.com

VERSATILE SINGLE-USE SANITARY FITTINGS INNOVATION SINGLE-USE SANITARY FITTINGS
LEADERSHIP SINGLE-USE SANITARY FITTINGS
TRUSTED SINGLE-USE SANITARY FITTINGS
RELIABLE SINGLE-USE SANITARY FITTINGS
ADVANCE FITTINGS

**SPECIALISTS
IN FLUID
CONNECTIONS**

Call, or email, today to find out how our Specialists in Fluid Connections can help you find a fitting for your application!

THE SF SERIES

SANITARY FITTINGS

*Economical Single-Use
Animal-Free Materials
High Flow Capacity
Mini & Maxi Sizes*



SFMN655

Mini Flange to 600 Series Barb
1/4" (6.4 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMN670

Mini Flange to 600 Series Barb
3/8" (9.5 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMN680

Mini Flange to 600 Series Barb
1/2" (12.7 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMN6100

Mini Flange to 600 Series Barb
3/4" (19.0 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMX655

Maxi Flange to 600 Series Barb
1/4" (6.4 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMX670

Maxi Flange to 600 Series Barb
3/8" (9.5 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMX680

Maxi Flange to 600 Series Barb
1/2" (12.7 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMX6100

Maxi Flange to 600 Series Barb
3/4" (19.0 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMX6110

Maxi Flange to 600 Series Barb
1.0" (25.4 mm) Tubing ID
Animal-Free Polypropylene
Natural PVDF



SFMNC

Mini Cap
Animal-Free Polypropylene



SFMXC

Maxi Cap
Animal-Free Polypropylene



GMN450

Molded Gasket for Mini Flange
Platinum-Cured Silicone



GMX455

Molded Gasket for Maxi Flange
Platinum-Cured Silicone

**SPECIALISTS
IN FLUID
CONNECTIONS**



VALUE PLASTICS, INC.
SPECIALISTS IN FLUID CONNECTIONS™

HROMalytic +61(0)3 9762 2034 09
AUSTRALIAN Distributors **ECHnology** www.chromtech.net.au

