



Ethanol Analytical Reference Standards for Blood Alcohol Testing

Capillary GC Columns



Rtx®-BAC1 (fused silica)

ID	df (µm)	30-Meter
0.32mm	1.80	18003
0.53mm	3.00	18001

Rtx®-BAC2 (fused silica)

ID	df (µm)	30-Meter
0.32mm	1.20	18002
0.53mm	2.00	18000

MXT®-BAC1 (Silcosteel®-treated stainless)

ID	df (µm)	30-Meter
0.53mm	3.00	78001

MXT®-BAC2 (Silcosteel®-treated stainless)

ID	df (µm)	30-Meter
0.53mm	2.00	78000


- NIST-traceable ethanol calibration standards from 0.015g/dL to 0.4g/dL in water.
- Eight-component resolution control standard with NIST-certified ethanol concentration.
- Datapack and Certificate of Analysis provided with each standard.

Legal limits for blood alcohol in the United States are set at either 0.080 or 0.100g/dL. We offer high- and low-concentration standards so that laboratories can construct calibration curves that bracket either legal limit concentration. In addition, we offer 0.015g/dL and 0.025g/dL concentrations, which are important in “zero-tolerance” cases involving minors or convicted DUI offenders.

These standards are National Institute of Standards and Technology (NIST)-traceable and are shipped with a complete data pack that includes a Certificate of Analysis, copies of all raw material testing results, statistical QA results, analytical balance printout, and a lot sheet showing gravimetric weight of each analyte.

Use our Resolution Control Standard mix to verify the retention time for each compound normally detected in a blood alcohol test, and to verify that the compounds are resolved from and do not interfere with each other. Each compound in this resolution mix is 0.100g/dL. The ethanol concentration is an NIST-traceable certified value.

We also offer Rtx®-BAC1 and Rtx®-BAC2 columns, which can resolve a blood alcohol sample in less than 3 minutes. If you have any questions about blood alcohol methods, or these products, please contact our Technical Service Team at support@restekcorp.com or via phone at 800-356-1688 or 814-353-1300, ext. 4. To order call ext. 3.

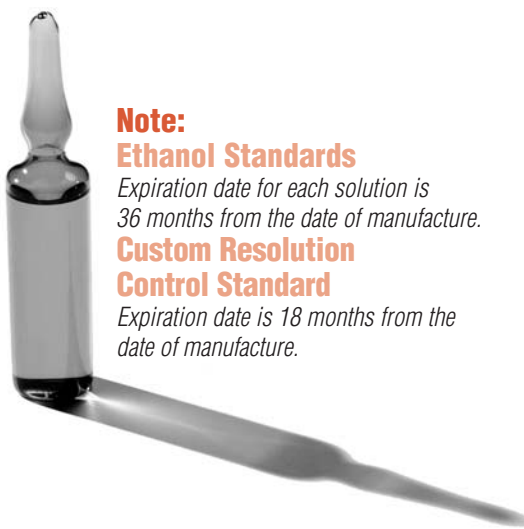

Note: Ethanol Standards

Expiration date for each solution is 36 months from the date of manufacture.

Custom Resolution Control Standard

Expiration date is 18 months from the date of manufacture.

Forensic Ethanol Solutions w/data pack

Forensic ethanol solutions w/data pack	5-pk. 1mL/ampul	10-pk. 1mL/ampul	ea. 5mL/ampul	ea. 20mL/ampul
0.015g/dL forensic ethanol solution	36232	36332	36240	36248
0.02g/dL forensic ethanol solution	36233	36333	36241	36249
0.025g/dL forensic ethanol solution	36234	36334	36242	36250
0.04g/dL forensic ethanol solution	36235	36335	36243	36251
0.05g/dL forensic ethanol solution	36257	36259	36258	36260
0.08g/dL forensic ethanol solution	36262	36264	36263	36265
0.1g/dL forensic ethanol solution	36236	36336	36244	36252
0.15g/dL forensic ethanol solution	36237	36337	36245	36253
0.2g/dL forensic ethanol solution	36238	36338	36246	36254
0.3g/dL forensic ethanol solution	36239	36339	36247	36255
0.4g/dL forensic ethanol solution	36266	36268	36267	36269

Blood Alcohol Mix Resolution Control Standard (8 components)

acetaldehyde	ethanol (NIST certified value)	methanol
acetone	ethyl acetate	methyl ethyl ketone
acetonitrile	isopropanol	

0.100g/dL each in water; 1mL/ampul

Each w/data pack
36256

RESTEK 800-356-1688
814-353-1300
www.restekcorp.com

HRMalytic www.chromtech.net.au sales@chromtech.net.au
ABN 14 643 445 058 PTY LTD Tel: (03) 9762 2034 Fax: +61 3 9761 1169
Australian Distributors: **TECH**nology . . . for ALL Your chromatography supplies!