


Thermal Desorption Unit (TDU) Tubes

Variety of sorbents to collect a wide range of VOCs



HROMalytic Chromatography
Australian Distributors **ECH**nology Products '08
www.chromtech.net.au E-mail : info@chromtech.net.au Tel : +61 3 9762 2034 Fax : +61 3 9761 1169

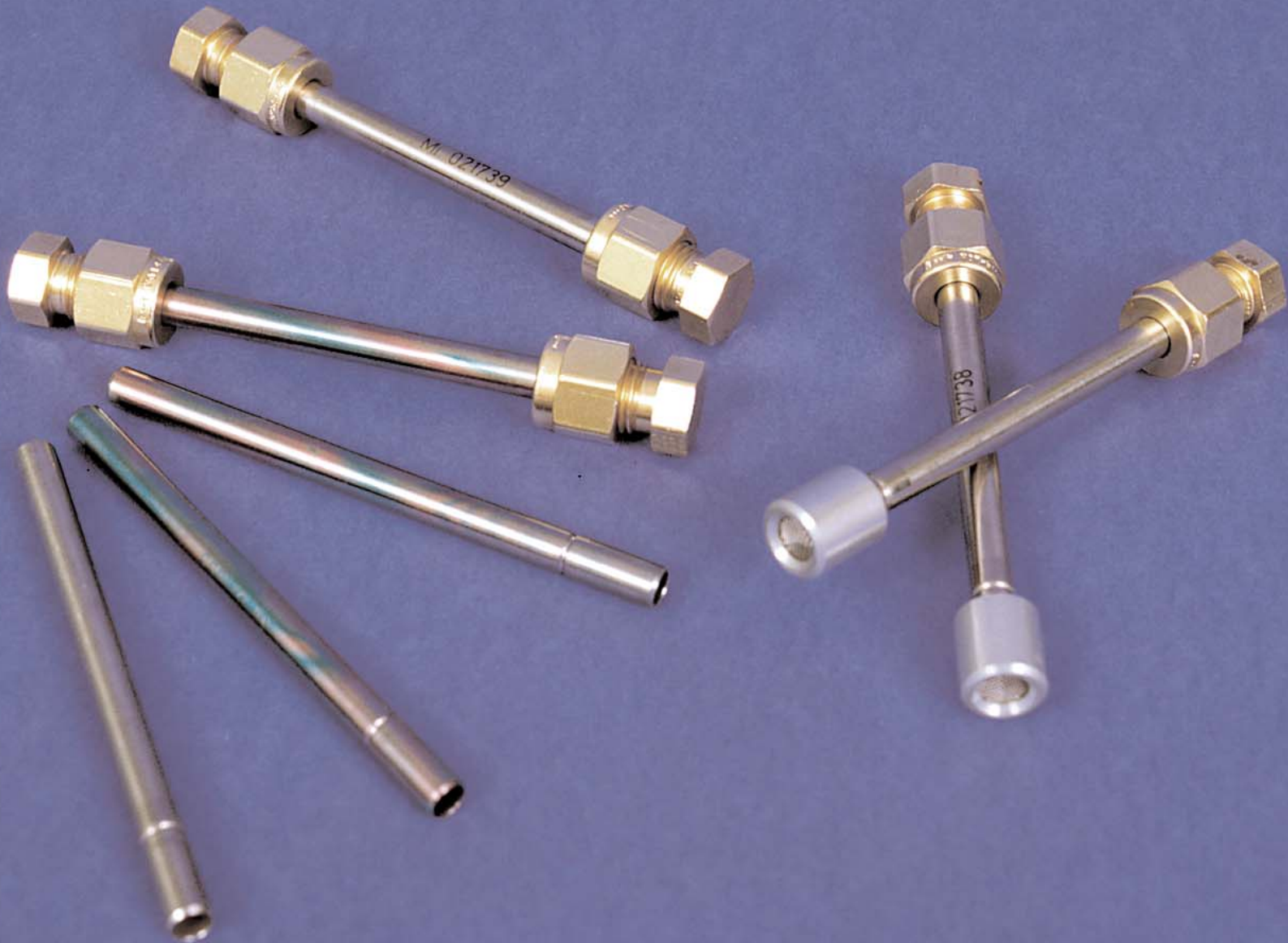

Chromatography Products

www.restek.com

Thermal Desorption Unit (TDU) Tubes

- Variety of sorbents to collect a wide range of VOCs.
- Use glass tubes for maximum inertness in active sampling.
- Choose stainless steel tubes for either active or passive sampling. No sampling pump necessary for passive sampling with diffusion caps!
- Individually etched with unique serial number for convenient sample identification.
- Available unconditioned or preconditioned and ready to sample. Tubes are reusable after thermal desorption.

High-quality thermal desorption tubes by Markes International are now available from Restek. These sorbent tubes are suitable for ppt to ppm concentrations of volatile organic compounds (VOCs) in ambient, indoor, and industrial hygiene environments. Available in both stainless steel and glass (for thermally labile VOCs), they fit Markes ULTRA-UNITY, PerkinElmer, and Shimadzu thermal desorbers. Packed tubes come with a report detailing the total mass of sorbent in the tube; conditioned tubes also include a blank chromatogram.



method applications

Method	Application
US EPA	TO-17
ASTM	D-6196
NIOSH	2549
DIN EN ISO	16017

Specifications

Dimensions: 1/4" OD x 3-1/2" long
 Low sampling rates:
 0.01-0.20 L/min.
 (<10L total volume)
 Long-term storage caps are supplied with conditioned tubes

Thermal Desorption Tube Sorbent	Applications	Unconditioned		Conditioned & Capped	
		Stainless Steel	Glass	Stainless Steel	Glass
Tenax TA	Vapor phase organics from C6/7 to C26				
Graphitized Carbon	Vapor phase organics from C5/6 to C14				
Tenax GR/Carbopack B	Vapor phase organics from <i>n</i> -C5/6 to <i>n</i> -C20 (EPA Methods TO-14/TO-15/TO-17)				
Carbopack B/Carbosieve SIII	Vapor phase organics from <i>n</i> -C2/3 to <i>n</i> -C12/14 (EPA Methods TO-14/TO-15/TO-17)				
Tenax TA/Graphitized Carbon/Carboxen 1000	Vapor phase organics from C2/3 to C20				
Carbopack C/Carbopack B/Carbosieve SIII	Vapor phase organics from <i>n</i> -C2/3 to <i>n</i> -C16/20 (EPA Methods TO-14/TO-15/TO-17)				

Description	qty.	Unconditioned		Conditioned & Capped	
		Stainless Steel cat.#	Glass cat.#	Stainless Steel cat.#	Glass cat.#
TDU Tubes, Tenax TA	10-pk.	24056	24062	24080	24086
TDU Tubes, Graphitized Carbon	10-pk.	24057	24063	24081	24087
TDU Tubes, Tenax GR/Carbopack B	10-pk.	24058	24064	24082	24088
TDU Tubes, Carbopack B/Carbosieve SIII	10-pk.	24059	24065	24083	24089
TDU Tubes, Tenax TA/Graphitized Carbon/Carboxen 1000	10-pk.	24060	24066	24084	24090
TDU Tubes, Carbopack C/Carbopack B/Carbosieve SIII	10-pk.	24061	24067	24085	24091



Stainless Steel, Conditioned and Capped

Thermal Desorption Unit Tubes, Empty

- Empty tubes for direct desorption of VOCs in liquids, solids, or pastes.
- Stainless steel: front sorbent retaining gauze fitted, rear gauze and gauze retaining spring supplied.
- Glass: with glass frit positioned 15mm from sampling end.

Description	qty.	Stainless Steel	Glass
		cat.#	cat.#
TDU Tubes, Empty	10-pk.	24054	24055



Glass, Unconditioned

Thermal Desorption Unit Tubes, Calibration

Description	qty.	Stainless Steel	Glass
		cat.#	cat.#
TDU Tubes, Calibration, Tenax TA 1cm Bed	10-pk.	24075	24076

Description	qty.	cat.	
		Stainless Steel	Glass
Calibration Solution Loading Rig	ea.	24077	
Calibration Solution Loading Rig 9.5mm Replacement Septa	10-pk.	24078	
Certified Reference Standard, 100ng BTX on Tenax TA	10-pk.	24079	



Stainless Steel, Unconditioned

Thermal Desorption Unit Tubes, Accessories

Description	Benefits/Uses	qty.	cat.
1/4" Brass Cap and PTFE Ferrules	Use for long-term storage of blank/sampled tubes.	20-pk.	24068
1/4" PTFE Ferrules	Long-term storage caps.	20-pk.	24069
CapLok Tool	Use for tightening long-term storage caps.	ea.	24070
Pen Clip		10-pk.	24071
TubeMate Tool	Assists with tube packing.	ea.	24072
1/4" Stainless Steel Union and PTFE Ferrules	Use for connecting tubes in series.	10-pk.	24073
Diffusion Caps	Required for diffusive sampling with stainless steel tubes.	10-pk.	24074



CapLok Tool



Diffusion Caps

Restek Trademarks: Carbopack Carbosieve, Restek logo.



www.chromtech.net.au E-mail : info@chromtech.net.au Tel : +61 3 9762 2034 Fax : +61 3 9761 1169