


Ana-Box™ Analog Sensor Interface-\$99

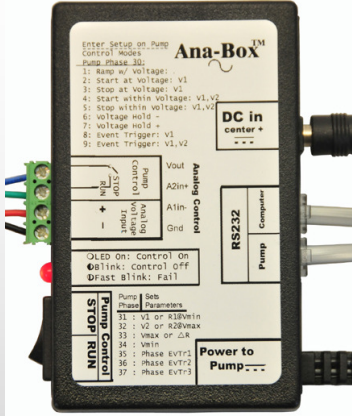

MARKET EXCLUSIVE


Interface your pressure sensors, PLCs, sensors, joysticks, etc. with the Ana-Box™ Analog Sensor Interface to your pump.

***So capable...
yet so incredibly simple to use!***


Sample Ana-Box™ Operational

- Over 10 modes of operation
- Vary the pumping rate in proportion to a variable voltage input
- Start or stop pump at set voltages
- Vary the pumping rate to maintain a pressure or pH level
- Attaches to any pump in the NE-1000 Family
- Operates stand alone. Does not need a computer to operate
- Trigger program events according to voltage
- Optionally outputs event data to a computer for data logging
- No complicated user interface.
- OEM customizing available


Ramping (gradient) with Voltage
Automatically adjusts the pumping rate proportionally to the voltage

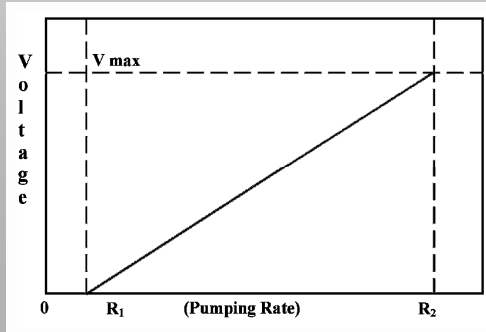


Start Pump at Voltage Trigger

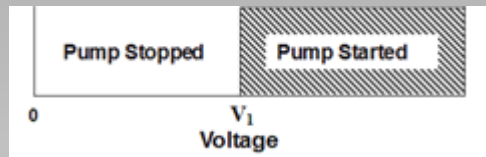


Stop Pump at Voltage Trigger

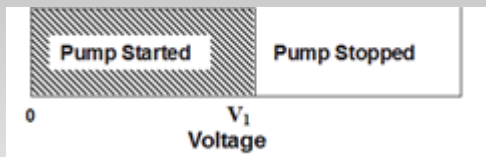


Pressure Controlled Infusion
Automatically varies the pumping rate to maintain a pressure by utilizing your pressure sensor with a voltage output.

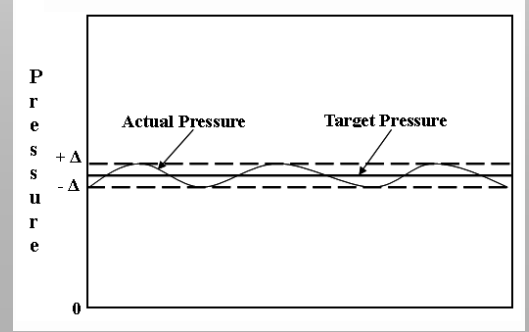

Start Pump within Voltage

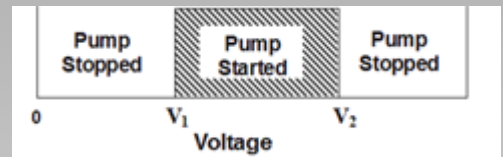

Stop Pump within Voltage

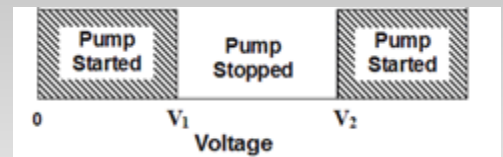

SyringePump.com
Clever Pumps, Priced Right!

138 Toledo St.
Farmingdale, NY 11735
Phone: (631) 249-1392
Fax: (707) 248-2089
info@syringepump.com



Australian Distributors
Importers & Manufacturers
www.chromtech.net.au

Website NEW : www.chromalytic.net.au E-mail : info@chromtech.net.au Tel: 03 9762 2034 . . . in AUSTRALIA