


Series A-2

Removable needles: .028" x .005" x 2", bevel, open end on 25, 50, and 100 μ l
.029" x .012" x 2", bevel, open end on all others

All the advantages of the Series A syringe plus:

- ▶ **Push-button valve** – allows instantaneous injection
- ▶ **Smaller volumes** – great for small volatile samples

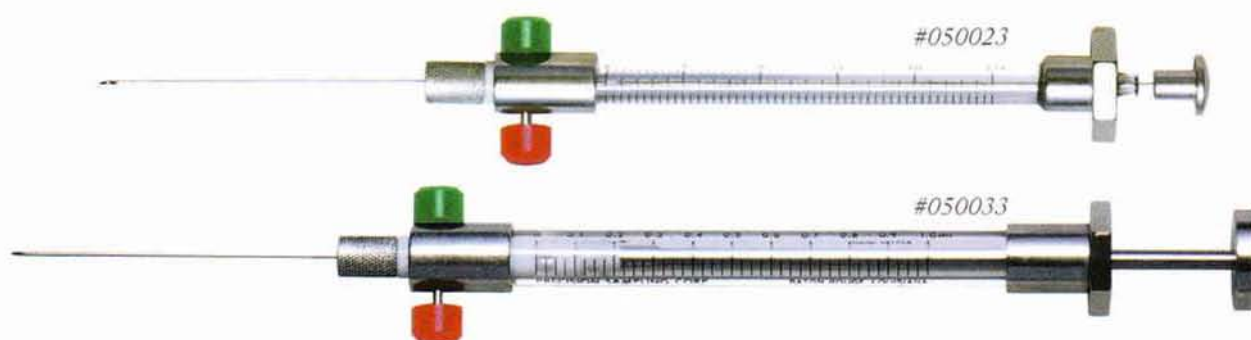

The A-2 features a push-button valve, which allows sample storage to 250 psi in syringes as small as 25 μ l. This is useful with small liquid samples containing low-boiling components which would be lost through evaporation using ordinary syringes.

The positive rear stop (in 250 μ l and larger sizes) prevents plunger blowout at elevated pressures, protecting against sample loss or operator injury. Like the Series A gas syringe, the A-2 has all the standard Pressure-Lok features such as a PTFE plunger tip, PTFE-sealed needle, and ultrasmooth bore. Replacement components (page 12) are available for easy repair.

	25 μ l	50 μ l	100 μ l	250 μ l	500 μ l
Standard:	050023	050024	050025	050031	050032
Luer:	050043	050044	050045	050051	050052

	1 ml	2 ml	5 ml	10 ml
Standard:	050033	050034	050035	050036
Luer:	050053	050054	050055	050056

Replacement needles for standard syringes (Luer needles on page 11)			
Size:	.028" x .005" x 2"	.029" x .012" x 2"	.029" x .012" x 2"
Type:	Bevel, open end	Bevel, open end	Side port, taper
Pkg. of 3:	943050	943051	943052

Gas Syringes

Series C

- ▶ **Optional Chaney adapter** – much lower dead volume than luer-style connectors
- ▶ **Perfect syringe for everyday use** – all our standard features, but without a valve for pressurized sample storage



	100 µl	250 µl	500 µl	1 ml	2 ml	5 ml	10 ml
Standard:	030025	030031	030032	030033	030034	030035	030036
Std. w/Chaney adapter:	101430	101440	101448	101455	101463	101470	101478
Luer:	030045	030051	030052	030053	030054	030055	030056

Replacement needles (Luer needles on page 11)

	Standard		Chaney adapter
Size:	.029" x .005" x 2"	.029" x .012" x 2"	.018" x .008" x .75"
Type:	Bevel, open end	Bevel, open end	Side port, taper
Pkg. of 3:	943050	943051	101935

Add sample storage capability to any syringe with a Mininert syringe valve