

Series D

Removable needle: .029" x .012" x 2", bevel, open end

An economically priced high precision gas sampling syringe which offers more features than conventional gas-tight designs: leak-tight, Pressure-Lok plunger tip; no dead volume barrel design; and special low dead volume needle connector for accurate gas samples.



	1 ml	2 ml	Replacement needles (Luer needles on page 11)	
Standard:	040033	040034	Size: .029" x .012" x 2"	.029" x .012" x 2"
Luer:	040053	040054	Type: Bevel, open end	Side port, taper
			Pkg. of 3: 943051	943052


Mininert Syringe Valves

For Series C and D syringes or conventional luer-tip syringes

- ▶ **Add sample storage capability** – mount a Mininert syringe valve on a standard syringe
- ▶ **Push button operation** – green to open, red to close

For a Series C or D syringe

These convenient add-on valves allow our Series C and D syringes to store samples at up to 250 psi. The valve body is all PTFE, with a stainless steel stem.

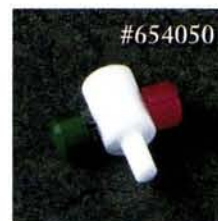

For a Luer-tip syringe

This model is identical to the valve above, except that it is designed to fit Luer-tip syringes from any manufacturer and accept traditional Luer needles.



For C or D syringe	For luer-tip syringe
654050	654051