

# ROTO-SHAKE GENIE®

Rock, Rotate, or 3-D Orbit Any Container!



More multi-task products from the makers of the Vortex-Genie®!

- Combines rocking and rotating in one compact and versatile unit
- Rocking tray accommodates large open containers
- Magnetic strips secure sealed bags for completely uniform soaking and washing
- Optional 3-D orbital attachment transforms the unit into an orbital shaker
- Easily cleaned or decontaminated
- **MAGNETIC** stainless steel platform may be easily set to rock or rotate almost any container
- A variety of attachments are available to provide maximum application versatility
- Small footprint and economical price complement any type of laboratory
- May be used in cold rooms and incubators between 0-38°C

U.S. Pat. 5,918,979

## Ordering Information

Catalog No.	Description
SI-1100	Roto-Shake Genie, 120V, 60Hz
SI-1101	Roto-Shake Genie, 230V, 50Hz, No Plug
SI-1102	Roto-Shake Genie, 230V, 50Hz, European Plug
SI-1103	Roto-Shake Genie, 230V, 50Hz, British Plug
SI-1104	Roto-Shake Genie, 230V, 50Hz, Swiss Plug
SI-1105	Roto-Shake Genie, 100V, 50/60Hz

## Specifications

Supplied	1 Clip Plate each for 10-13 mm, 15-17mm, 28-30mm Tubes, 4 Bag Mounting Strips, Non-Skid Tray
Rotating	5-35 RPM
Rocking	10-70 cycles per minute
Rocking Angle	±10 degrees from horizontal
Platform Capacity	4.5 Kg (10 lb)
Weight	11.8 Kg (26 lb)
Dimensions (D x W x H)	252 x 508 x 267mm (9.9 x 20 x 10.5in)
Depth with Rocker Tray	305mm (12in)
Environmental	Installation Category 2, 0° - 38°C (32° - 100°F), 95% max. relative humidity




ETL Listed Conforms to UL Std. 3101-1  
ETL Listed Certified to CAN/CSA C22.2



www.chromtech.net.au E-mail : info@chromtech.net.au Tel : +61 3 9762 2034 Fax : +61 3 9761 1169

# Roto-Shake Genie® Accessories

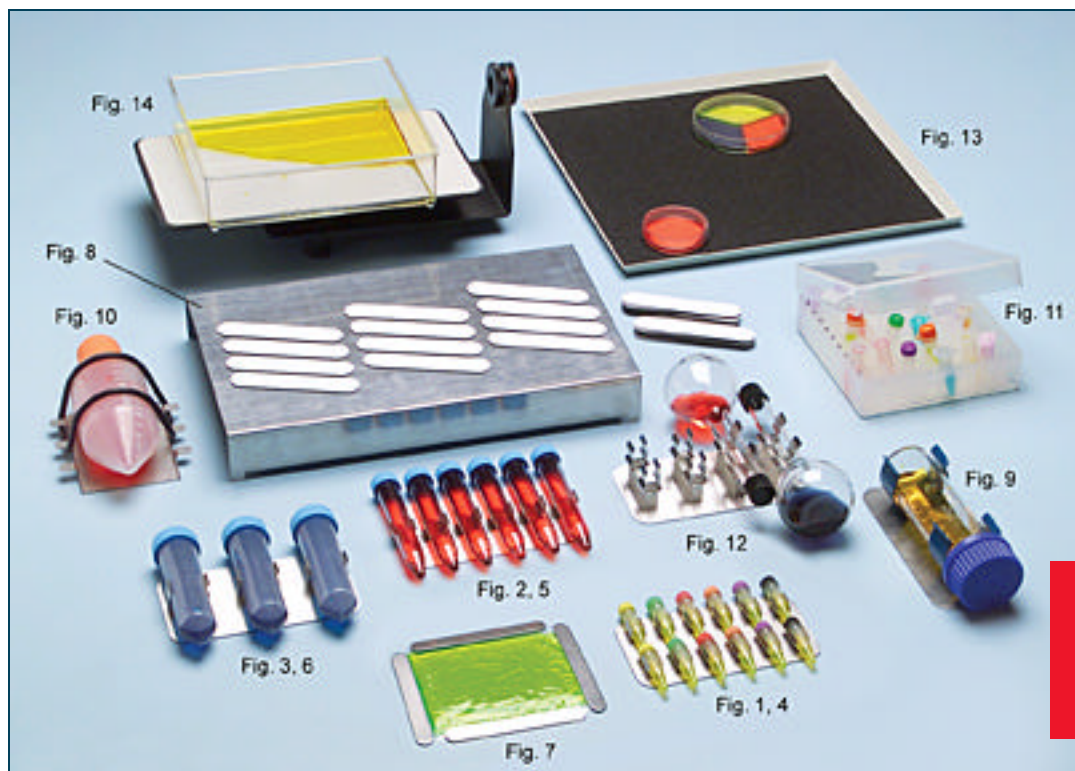

**Accessories  
attach  
magnetically  
to platform!**

**Accessories  
supplied  
with the  
Roto-Shake Genie**

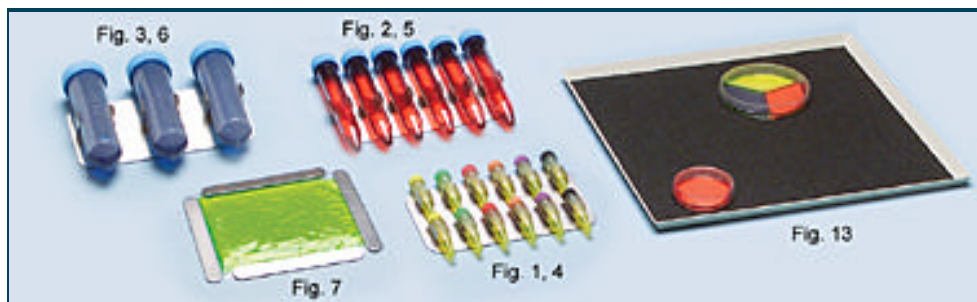


Fig.	Catalog No.	Description
1	SI-1123	1 Clip Plate for 12 each 10-13mm tubes
2	SI-1124	1 Clip Plate for 6 each 15-17mm tubes
3	SI-1125	1 Clip Plate for 3 each 28-30mm tubes
4	SI-1120	6 Clip Plates for 12 each 10-13mm tubes
5	SI-1121	6 Clip Plates for 6 each 15-17mm tubes
6	SI-1122	6 Clip Plates for 3 each 28-30mm tubes
7	SI-1126	Bag Mounting Strips (pack of 16)
8	SI-1127	Expansion Kit - two stainless steel trays, 24 bag mounting strips
9	SI-1130	2 Clip Plates for Hybridization Tubes
10	SI-1134	2 Universal Clip Plates
11	SI-1135	1 Magnetic Covered 100 Microtube Box
12	SI-1136	1 Clip Plate for 100ml Volumetric Flasks
13	OA-1100-070	Non-Skid Tray
14	SI-1150	3-D Orbital Shaker Attachment



*Plasticware, glassware and other consumable products shown are not included.*



www.chromtech.net.au E-mail : info@chromtech.net.au Tel : +61 3 9762 2034 Fax : +61 3 9761 1169