

FID Replacement Jets

Standard Version

- Engineered with a fluted tip to guide the capillary column into the jet.
- Threads specially coated for easy installation and removal.
- Special processing ensures the highest degree of cleanliness.

High-Performance Version

- Similar to the standard version, but Siltek® treated.
- Extremely inert, for use with active compounds.

Capillary Adaptable FID Replacement Jet for Agilent 5890/6890/6850 GCs

0.011-Inch ID Tip	Similar to Agilent part #	qty.	cat.#	qty.	cat.#
Standard, 0.011-Inch ID Tip	19244-80560	ea.	20670	3-pk.	20671
High-Performance Siltek® Treated, 0.011-Inch ID Tip	19244-80560	ea.	20672	3-pk.	20673

Capillary Dedicated FID Replacement Jet for Agilent 6890/6850/7890 GCs

0.011-Inch ID Tip	Similar to Agilent part #	qty.	cat.#	qty.	cat.#
Standard, 0.011-Inch ID Tip	G1531-80560	ea.	21621	3-pk.	21682
High-Performance Siltek® Treated, 0.011-Inch ID Tip	G1531-80560	ea.	21620	3-pk.	21683

Packed Column FID Replacement Jets for Agilent 5890/6890/6850 GCs

0.018-Inch ID Tip	Similar to Agilent part #	qty.	cat.#	qty.	cat.#
Standard, 0.018-Inch ID Tip	18710-20119	ea.	21694	3-pk.	21695
High-Performance Siltek® Treated, 0.018-Inch ID Tip	18710-20119	ea.	21696	3-pk.	21697

0.030-Inch ID Tip	Similar to Agilent part #	qty.	cat.#	qty.	cat.#
Standard, 0.030-Inch ID Tip	18789-80070	ea.	21688	3-pk.	21689
High-Performance Siltek® Treated, 0.030-Inch ID Tip	18789-80070	ea.	21686	3-pk.	21687




tech tip

Which FID Jet Should I Use?

There are two FID jet configurations for Agilent GCs. The longer "adaptable" jet fits both 5890 and 6890 GCs, and can be used with capillary or packed columns. The shorter "dedicated" jet is for the FID in the 6890 GC that is designed only for use with capillary columns.

FID Jet Removal Tool for Agilent 5890/6890/6850/7890 FIDs

- Securely grips jet in socket for easy removal or installation.
- Unique, ergonomic handle—easy to hold.

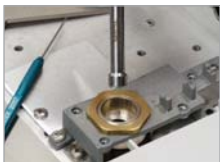

Slip tool over FID jet...



...loosen jet...



...and remove.

Description	qty.	cat.#
FID Jet Removal Tool for Agilent 5890/6890/6850/7890 FIDs	ea.	22328

restek
innovation!



Torx® Screwdriver Set

- Set includes TR-10, TR-15, and TR-20.
- Ideal for performing routine maintenance on Agilent 6890 and 7890 GCs.

Description	qty.	cat.#
Torx® Screwdriver Set	set	23034

