

Supplies for Agilent Instruments: Detector Replacement Parts



Detector Plug Nut

Use Restek's detector plug nut to cap off a detector for thermal cleaning, to check detector or make-up gas flow rates, or to prevent hydrogen from accidentally diffusing into the oven from an unused detector base.

Description	Similar to Agilent part #	qty.	cat. #
Detector Plug Nut	5020-8294	2-pk.	21883



FID/NPD Adaptor Fitting

- Easy to use, sturdy, compact stainless steel fitting.
- 1/16-inch nut uses standard graphite or Vespel®/graphite ferrules.
- Wrench pad won't turn when installing a capillary column.
- Includes 1/4- and 1/16-inch stainless steel nuts, and 1/4-inch Vespel® and 0.4mm ID graphite ferrules.


Description	qty.	cat. #
FID/NPD Adaptor Fitting	kit	20884



FID/NPD Capillary Adaptor Fitting for Agilent 5890/6890/6850 GCs

- High-quality stainless steel construction.
- Meets or exceeds manufacturer's performance.

Description	Similar to Agilent part #	qty.	cat. #
FID/NPD Capillary Adaptor Fitting	19244-80610	ea.	22077

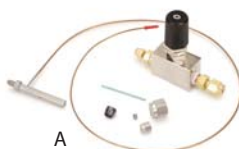

ECD/FID Dual-Purpose Make-Up Gas Kit for Agilent 5890 GCs

Incorporates a better make-up gas flow profile and a straight section of inert Siltek® tubing as the capillary guide. To add make-up gas to an FID, remove the tubing guide and insert an FID jet tail. A direct replacement fitting with a special end connector is available for GCs factory-equipped with make-up gas. A complete kit is available for plumbing GCs that are not factory-equipped with make-up gas.

Description	Similar to Agilent part #	qty.	cat. #
A) ECD/FID Dual-Purpose Make-Up Gas Kit Includes: make-up gas fitting, Siltek® treated guide, 1/4" nut and ferrule, 1/16" nut, 0.4mm ID graphite ferrule, and fine metering valve		kit	21300
B) Replacement Make-up Gas Fitting Kit Includes: make-up gas fitting, 1/4" nut and ferrule, 1/16" nut, 0.4mm ID graphite ferrule		kit	21324
C) ECD/FID Replacement Fitting Kit with Flow Manifold Connection Includes: replacement fitting, 1/4" nut, Vespel®/graphite ferrule, 1/16" nut, 0.4mm ID graphite ferrule, Siltek® treated guide		kit	21301
Replacement ECD Siltek® Metal Liner		2-pk.	21302
Replacement ECD Fused Silica Liner	19233-20625	ea.	27192
Replacement Micro ECD Liner	G2397-20540	ea.	27190