

Resprep™ SPE Manifolds, Vacuum/Pressure Unit




Resprep™ 12-Port SPE Manifold

Resprep™ 12- or 24-Port SPE Manifolds

- Use with any standard male luer-end SPE cartridges.
- Inert, Teflon® sample guides reduce cross-contamination and carryover.
- Flexible sample collection rack will accommodate a variety of receiving vessels.
- Quick vacuum-release valve for better system control.
- Individual vacuum control for each cartridge—improves reproducibility.

SAMPLE PREPARATION

Resprep™ Manifold Replacement Parts

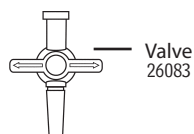

Valve
26083




Valve
Attachment
26130

Manifold Lid



Needle
Attachment
26131

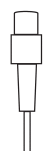

Sample
Guide
26084



Resprep™ 12-Port SPE Manifold

Replacement Manifold Lid, 12-port

Collection Rack, 12-port

Description

Description	qty.	cat.#
Complete Resprep™ 12-Port SPE Manifold (includes glass basin with built-in vacuum regulator, polypropylene lid with 12 individual control valves, 12-position collection rack, 12 Teflon® sample guides, and waste container)	kit	26077
Complete Resprep™ 24-Port SPE Manifold (includes glass basin with built-in vacuum regulator, polypropylene lid with 24 individual control valves, 24-position collection rack, and 24 Teflon® sample guides)	kit	26080

Resprep™ Manifold Replacement Parts and Accessories

Description	qty.	cat.#
Replacement Waste Container, 12-port new!	ea.	24014
Collection Rack		
Collection Rack, 12-port	ea.	26079
Collection Rack, 24-port	ea.	26082
Manifold Lid		
Replacement Manifold Lid, 12-port (sample guides not included)	ea.	26078
Replacement Manifold Lid, 24-port (sample guides not included)	ea.	26081
Manifold Lid Replacement Parts		
Valves, 12 or 24-port	2-pk.	26083
Valve Attachment	48-pk.	26130
Needle Attachment	48-pk.	26131
Sample Guides, 12 or 24-port	12-pk.	26084

*Not included in 24-port manifold (cat.# 26080 or 26082).



Oil-Free Vacuum/Pressure Station for 12-Port Resprep™ SPE Systems

Allows switching from pressure to vacuum in a matter of seconds. Quiet, oil-free unit won't contaminate the extraction system or your samples. Provides 20" Hg (68kPa) vacuum or 18psi (124kPa) pressure.

Description	qty.	cat.#
Oil-Free Vacuum/Pressure Station (115V)	ea.	24002
Oil-Free Vacuum/Pressure Station (220V)	ea.	24003
Vacuum Tubing (10ft./3m, 1/4" ID)	ea.	24016

Not recommended for use with 24-port manifold; not available for 220VAC power outlets.