



Tetracycline Mix on Allure® Biphenyl

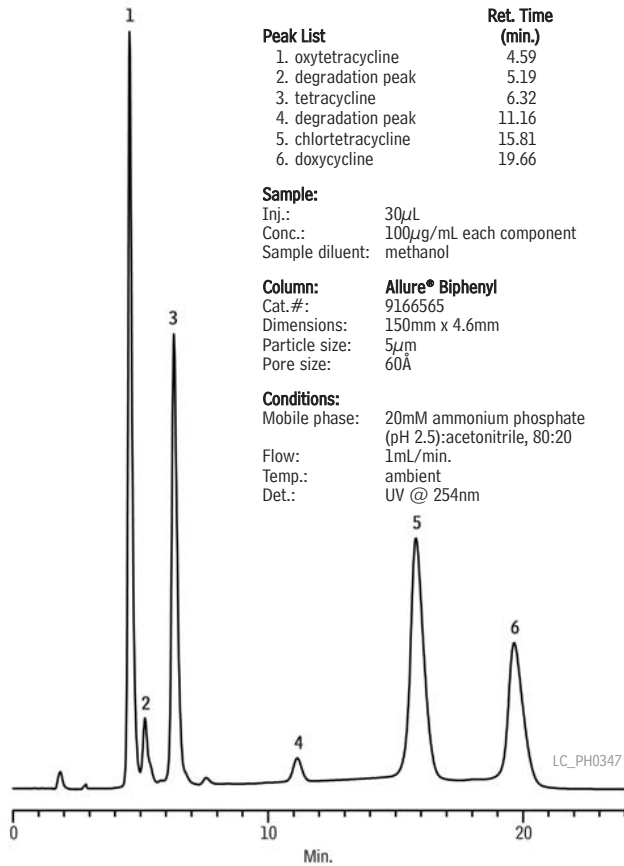

Oxytetracycline and Tetracycline on Allure® Biphenyl

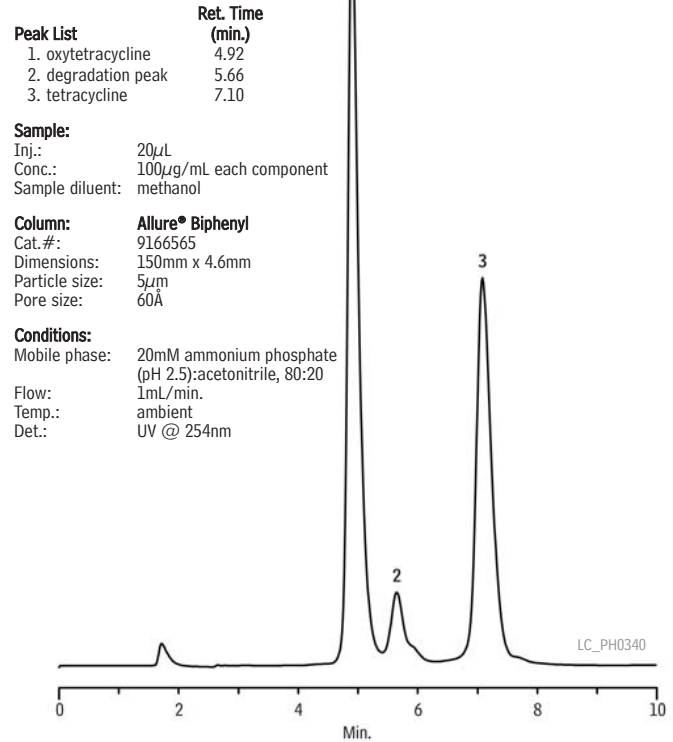

Amitriptyline on Ultra IBD, pH 7

new!


Peak List

1. amitriptyline

Sample:
Inj.: 10µL
Conc.: ~100µg/mL
Sample diluent: mobile phase

Column: Ultra IBD
Cat.#: 9175565
Dimensions: 150mm x 4.6mm
Particle size: 5µm
Pore size: 100Å

Conditions:
Mobile phase: 20mM potassium phosphate (pH 7):methanol, 10:90
Flow: 1.0mL/min.
Temp.: ambient
Det.: UV @ 254nm

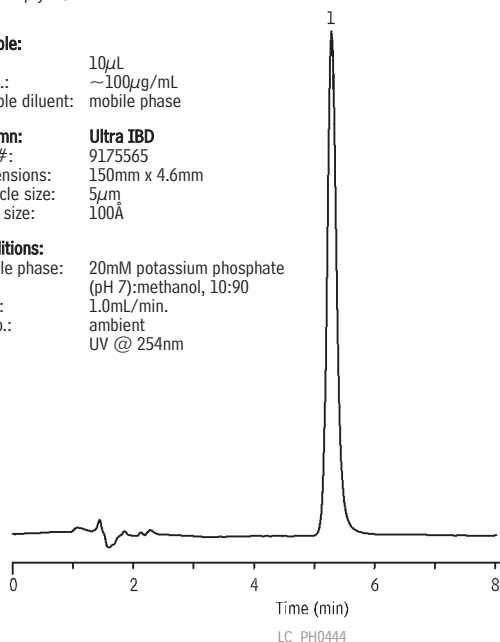

Searching for a chromatogram?
www.restek.com