



Aroclor PCBs (Quick Screening)
Rtx®-5


Aroclor 1221 Mix (cat.# 32007)

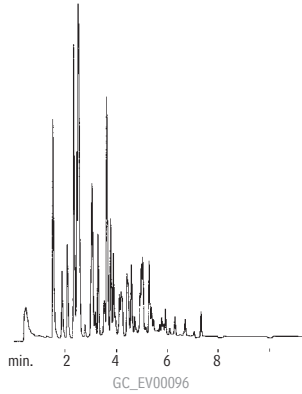


Aroclor 1232 Mix (cat.# 32008)

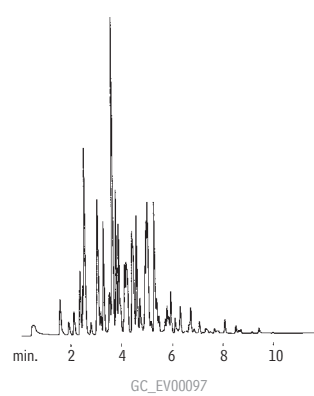


Aroclor 1242 Mix (cat.# 32009)

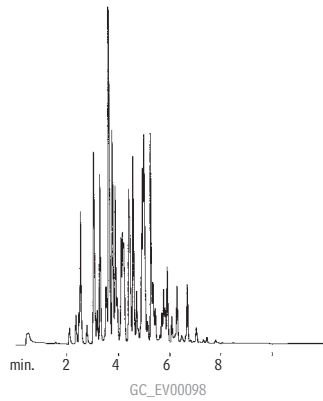


Aroclor 1248 Mix (cat.# 32010)

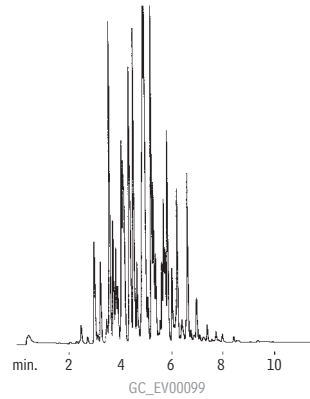

Aroclor 1254 Mix (cat.# 32011)

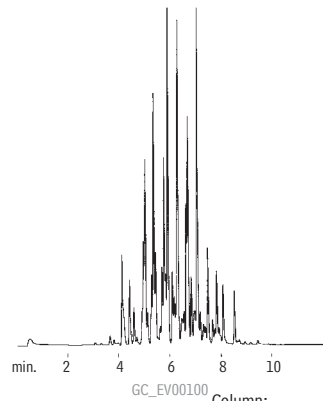

Aroclor 1260 Mix (cat.# 32012)

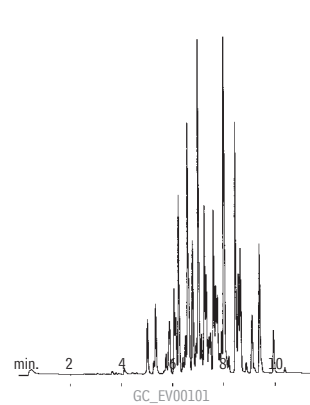

Column: Rtx®-5, 15m, 0.53mm ID, 0.50µm (cat.# 10237)
 Inj.: 1.0µL direct injection
 Conc.: 50ppm
 Oven temp.: 150°C to 300°C @ 12°C/min. (hold 5 min.)
 Inj./det. temp.: 250°C/310°C
 Carrier gas: helium
 Linear velocity: 30cm/sec. set @ 150°C
 ECD sensitivity: 16kHz full scale