


## Sugars (Alditol Acetates)

Rt™-2330

1. glyceraldehyde
2. deoxyribitol
3. rhamnitol
4. fucitol
5. ribitol
6. arabinitol
7. xylitol
8. deoxyglucitol
9. mannitol
10. galactitol
11. glucitol
12. inositol
13. glucoheptitol
14. *n*-acetyl galactose amine
15. *n*-acetyl glucose amine
16. 2-keto-3-deoxyoctanate

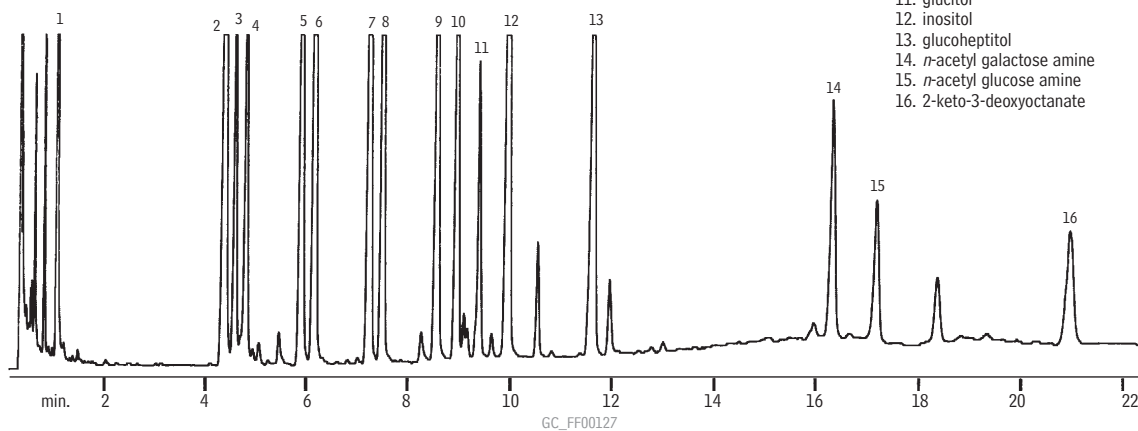

Column: Rt™-2330, 15m, 0.32mm ID, 0.20 $\mu$ m (cat.# 10721)  
 Inj.: 0.6 $\mu$ L split injection  
 Oven temp.: 175°C (hold 2 min.) to 240°C @ 8°C/min. (hold 1 min.) to 265°C @ 8°C/min. (hold 12 min.)  
 Inj./det. temp.: 275°C  
 Carrier gas: helium  
 Linear velocity: 80cm/sec. (flow rate: 10cc/min.)  
 FID sensitivity: 2 x 10<sup>-11</sup> AFS  
 Split ratio: 20:1

# Get More!

Food, Flavor & Fragrance  
 Related Articles Online

[www.restek.com/FFF](http://www.restek.com/FFF)