


Sugars (Alditol Acetates) Rtx®-225



Column: Rtx®-225, 15m, 0.25mm ID, 0.25 μ m (cat.# 14020)
 Inj.: 0.5 μ L split injection
 Oven temp.: 190°C (hold 5 min.) to 250°C @ 8°C/min. (hold 5 min.)
 Inj./det. temp.: 260°C
 Carrier gas: hydrogen
 Linear velocity: 42cm/sec. set @ 40°C
 FID sensitivity: 16 x 10¹¹ AFS
 Split ratio: 50:1

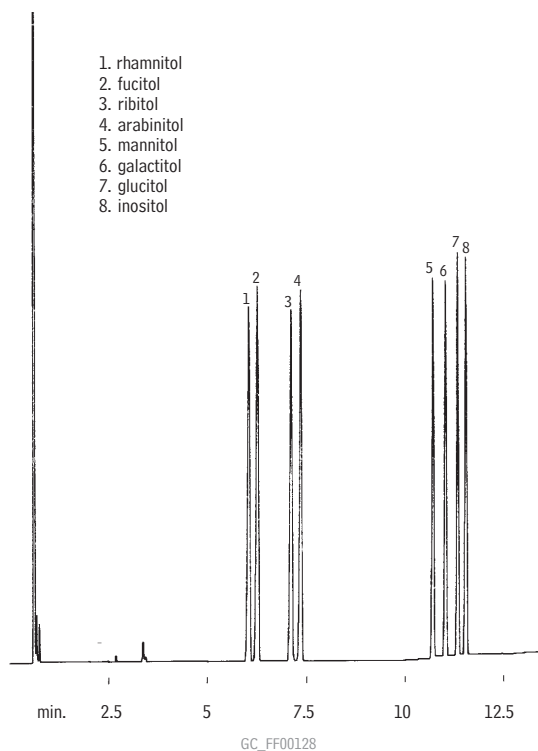

free literature

Foods, Flavors, and Fragrances Catalog!

Includes important analysis tips, and chromatograms for analysis of fats and oils, carbohydrates, vitamins, amino acids, organic acids, preservatives, flavors and fragrances, essential oils, and chiral separations. Retention time indices and complete product listings for all relevant GC and HPLC products also are included.
 lit. cat.# 59260A

Also available—Monitoring Volatile Compounds in Food Contact Packaging, Using Purge and Trap GC/MS and an Rtx®-5MS Capillary Column

Applications Note
 lit. cat.# 59348

Searching for free technical literature?

www.restek.com