

Trident Direct Guard Cartridge System

Easy to Use, Low Dead Volume—The Ultimate Combination of Convenience and Column Protection

Unlike “one size fits all” guard systems, the Trident Direct system gives you the power to select the right level of protection for your analysis. The system offers three levels of protection and guard cartridges in four dimensions, with a variety of bonded phases to match your analytical column. The economical, leak-free cartridge design provides an unprecedented combination of convenience, economy, and reliability. The foundation of the Trident Direct system is a reusable direct connect holder that easily attaches to any HPLC column using CPI- or Waters-style end fittings.* The system is available in configurations to match different protection level needs: in-line filter, in-line filter with holder for 10mm guard cartridge, and in-line filter with holder for 20mm guard cartridge. The guard cartridges are available in 2.1 and 4.0mm ID and are interchangeable within the appropriate length holder.



25082
Protection against particulate matter.



25084
Protection against particulate matter and moderate protection against irreversibly adsorbed compounds.



25086
Protection against particulate matter and maximum protection against irreversibly adsorbed compounds.

Description	qty.	cat.#	price
High-pressure filter	ea.	25082	
10mm guard cartridge holder without filter	ea.	25083	
10mm guard cartridge holder with filter	ea.	25084	
20mm guard cartridge holder without filter	ea.	25085	
20mm guard cartridge holder with filter	ea.	25086	
Connection tip for Waters-style end fittings	ea.	25088	
PEEK tip standard fittings	ea.	25087	
Replacement Cap Frits: 4mm 2.0µm	5-pk.	25022	
Replacement Cap Frits: 4mm 0.5µm	5-pk.	25023	
Replacement Cap Frits: 2mm 2.0µm	5-pk.	25057	
Replacement Cap Frits: 2mm 0.5µm	5-pk.	25990	

*The standard PEEK tip in Trident Direct systems is compatible with Parker, Upchurch Scientific, Valco, and other CPI-style fittings. To use Trident Direct systems with Waters-style end fittings, replace the tip with cat.# 25088.

Trident HPLC In-Line Guard Cartridge Holders

A Trident in-line guard cartridge holder can be used with almost any HPLC column by connecting it with a short piece of 1/16" tubing, appropriate nuts and ferrules, or finger-tight fittings. The system can be used with Restek columns, or with columns from other manufacturers. Holders are available for either 10mm or 20mm guard cartridges. Either size can be purchased with or without a prefilter, which provides added protection against the particles that can shorten the lifetime of the guard cartridge.



25021



25040



25061




25060

Description	qty.	cat.#	price
Holder for 10mm guard cartridge	ea.	25021	
Holder with filter for 10mm guard cartridge	ea.	25040	
Holder for 20mm guard cartridge	ea.	25061	
Holder with filter for 20mm guard cartridge	ea.	25060	
Replacement Cap Frits: 4mm 2.0µm**	5-pk.	25022	
Replacement Cap Frits: 4mm 0.5µm	5-pk.	25023	
Replacement Cap Frits: 2mm 2.0µm**	5-pk.	25057	
Replacement Cap Frits: 2mm 0.5µm	5-pk.	25990	

**Standard porosity.

Trident Direct 10mm guard cartridge holder with filter

Components

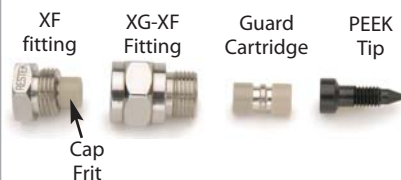

Assembled




Installed onto column



Cap frits

Trident In-Line 10mm guard cartridge holder with filter

Components

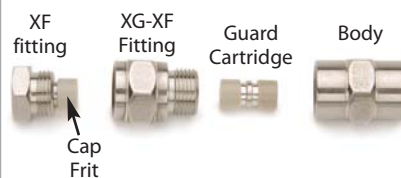

Assembled



Installed onto column

