




Brass



Stainless Steel

### Capillary Column Nuts

for Agilent 5890/6890/7890 GCs  
Available in brass or stainless steel.

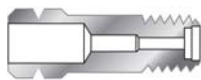

for use with "compact"  
Agilent-style ferrules



for use with standard ferrules

Description	Similar to Agilent part #	qty.	cat.#	price
<b>For use with "compact" Agilent-style ferrules.</b>				
Stainless Steel Capillary Column Nut	5181-8830	2-pk.	21884	
Brass Capillary Column Nut	5181-8830	2-pk.	21878	
<b>For use with standard 1/16"-type ferrules.*</b>				
Stainless Steel Capillary Column Nut	05921-21170	2-pk.	20883	
Brass Capillary Column Nut	05921-21170	2-pk.	21879	

\*Designed to fit a wider variety of 1/16" ferrules

### Hot Swap Capillary Column Nuts


- No more burned fingers!
- No more downtime waiting for injector parts to cool down.




22347

Never worry again about burned fingers or having to wait for the injector to cool down. The Hot Swap Capillary Column Nut allows you to change your capillary column while the injector temp is still hot.

Description	qty.	cat.#	price
<b>For use with "compact" Agilent-style ferrules.</b>			
Hot Swap Capillary Column Nut	ea.	22348	
<b>For use with standard 1/16"-type ferrules.</b>			
Hot Swap Capillary Column Nut	ea.	22347	

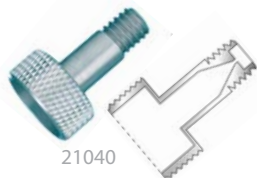

21040

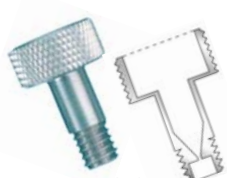

21041

### Finger-Tight Capillary Column Nuts

- Allow wrench-free column installations.
- Work with standard or compact (Agilent-style) ferrules.
- Made from high-quality stainless steel.

Description	Similar to Agilent part #	qty.	cat.#	price
<b>For use with "compact" Agilent-style ferrules.</b>				
Finger-Tight Capillary Column Nuts	5020-8293, 5020-8292	2-pk.	21040	
<b>For use with standard 1/16"-type ferrules.</b>				
Finger-Tight Capillary Column Nuts	5020-8293, 5020-8292	2-pk.	21041	



21312

### Finger-Tight Capillary Column Nuts

- Rapidly tighten without wrenches; avoid stripped threads.
- Either version can be used with 0.25, 0.32, or 0.53 mm ID columns.

Description	qty.	cat.#	price
<b>For use with "compact" Agilent-style ferrules.</b>			
Finger-Tight Capillary Column Nut	ea.	21311	
<b>For use with standard 1/16"-type ferrules.</b>			
Finger-Tight Capillary Column Nut	ea.	21312	



21034



21399

### Capillary Installation Gauge

- Seats graphite\* ferrules onto column for consistent installations.
- Prevents crushed column ends.
- Made from high-quality stainless steel.

Description	qty.	cat.#	price
Capillary Installation Gauge for Agilent-style fittings (compact ferrules)	ea.	21034	
Capillary Installation Gauge for 1/16" fittings (1/16" ferrules)	ea.	21399	

\*Recommended for use with graphite ferrules.

### Inlet Maintenance Kit for Agilent GCs

- Includes the most common consumable GC supplies and tools.
- All parts meet or exceed performance of instrument manufacturer's parts.
- Parts list makes reordering easy.



22181

#### Inlet kit includes:

- 0.4, 0.5, and 0.8 mm ID graphite ferrules.
- Viton® O-rings.
- Capillary nuts.
- Inlet seals.
- Reducing nut.
- Scoring wafer.
- 11 mm ThermoLITE® septa.
- 4.0 mm single gooseneck liner.
- 4.0 mm split liner with wool.
- Capillary column caps.
- 1/4" x 5/16" wrench.
- Septum puller.
- Installation gauge.
- Wire cleaning brush.
- Jet reamers/ferrule removers.
- Inlet liner removal tool.
- Septa nut removal tool.

Description	qty.	cat.#	price
Inlet Maintenance Kit	kit	22181	