

Torx® Screwdriver Set

- Set includes TR-10, TR-15, and TR-20.
- Ideal for performing routine maintenance on Agilent 6890 and 7890 GCs.



Description	qty.	cat.#	price
Torx Screwdriver Set	set	23034	



FID/Injector Cleaning Kit

The FID/Injector Cleaning Kit includes:

- Nylon tube brushes (1/8", 3/16", 1/4").
- Pipe cleaner.
- Stainless steel tube brushes (3/8", 3/16", 1/4").
- Stainless steel surface brush.
- Stainless steel jet reamers.
- Emery cloth.

Description	qty.	cat.#	price
FID/Injector Cleaning Kit	kit	20120	



Injector/Detector Plug Nuts

- Use to cap off an injector to isolate leaks.
- Use to cap off a detector for thermal cleaning.
- Use to check a detector or make-up gas flow rate.
- Use to cap off a detector and prevent hydrogen from accidentally diffusing into the oven from an unused detector base.



Description	Similar to Agilent part#	qty.	cat.#	price
Injector/Detector Plug Nuts	5020-8294	2-pk.	21883	

also **available**

Looking for more detector replacement parts?
See page 271.



FID/NPD Adaptor Fitting

- Easy to use, sturdy, compact stainless steel fitting.
- 1/16-inch nut uses standard graphite or Vespel®/graphite ferrules.
- Wrench pad won't turn when installing a capillary column.
- Includes 1/4- and 1/16-inch stainless steel nuts, and 1/4-inch Vespel® and 0.4mm ID graphite ferrules.

Description	qty.	cat.#	price
FID/NPD Adaptor Fitting	kit	20884	




FID/NPD Capillary Adaptor Fitting

for Agilent 5890/6890/6850 GCs

- High-quality stainless steel construction.
- Meets or exceeds original manufacturer's performance.

Description	Similar to Agilent part #	qty.	cat.#	price
FID/NPD Capillary Adaptor Fitting	19244-80610	ea.	22077	

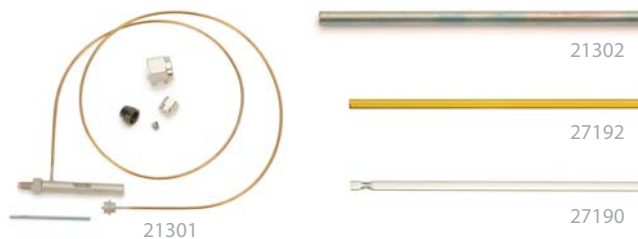

ECD/FID Dual-Purpose Make-Up Gas Kit

for Agilent 5890 GCs

Kit includes: replacement fitting, 1/4" nut, Vespel/graphite ferrule, 1/16" nut, 0.4 mm ID graphite ferrule, Siltek treated guide

Description	Similar to Agilent part #	qty.	cat.#	price
ECD/FID Replacement Fitting Kit with Flow Manifold Connection		kit	21301	
Replacement ECD Siltek Metal Guide		2-pk.	21302	
Replacement ECD Fused Silica Liner	19233-20625	ea.	27192	
		5-pk.	27193	
Replacement Micro ECD Liner	G2397-20540	ea.	27190	
		5-pk.	27191	