

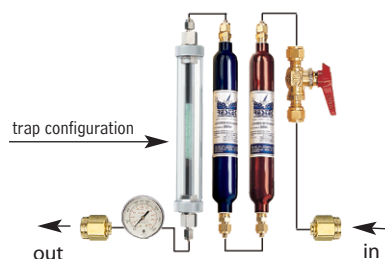
Restek Gas Management System

- Removes moisture, hydrocarbons, and oxygen from carrier gas, extending column lifetime.
- Produces high-purity carrier gas for most applications.
- Includes one each: moisture, hydrocarbon, and indicating oxygen trap.
- Replacing traps is safe and easy.
- Maximum flow: 1 liter/minute.



Dimensions:
12" x 14" x 3"
(30.5 x 35.6 x 7.6 cm)

Restek has put together a convenient unit providing gas purification all in one step. Complete with an indicating trap, your gas purification issues are handled in one central location.

**did you know?**

The Restek Gas Management System removes water vapor (to 10 ppb), hydrocarbons (to 0.1 ppm), and oxygen (to less than 0.1 ppm) with three traps housed in one unit.

Description	Fittings	qty.	cat.#	price
Restek Gas Management System	includes fittings for 1/8" and 1/4" gas line	kit	21999	
Replacement Traps		qty.	cat.#	price
High-Capacity Moisture Trap	1/8" Nickel-Plated Brass	ea.	21997	
Capillary-Grade Hydrocarbon Trap	1/8" Nickel-Plated Brass	ea.	21991	
Indicating Oxygen Trap	1/8" Nickel-Plated Brass	ea.	22010	

Restek Electronic Leak Detector

Don't let a small leak turn into a costly repair—protect your instrument and analytical column by using a Restek Leak Detector.

Backed by a 1-year warranty, the new Restek Leak Detector sets an industry standard for performance and affordability in hand-held leak detectors.

**Leak Detector Facts**

Detectable gases:	helium, nitrogen, argon, carbon dioxide, hydrogen
Battery:	Rechargeable Ni-MH internal battery pack (6 hours normal operation)
Operating Temperature Range:	32°-120°F (0°-48°C)
Humidity Range:	0-97%
Warranty:	one year
Certifications:	CE, Ex, Japan
Compliance:	WEEE, RoHS

Description	qty.	cat.#	price
Leak Detector with Hard-Sided Carrying Case and Universal Charger Set (US, UK, European, Australian)	ea.	22839	
Leak Detector Routine Maintenance Review**	ea.	22839-R	
Soft-Side Storage Case	ea.	22657	
Small Probe Adaptor	ea.	22658	

*Caution: The Restek Electronic Leak Detector is designed to detect trace amounts of hydrogen in a noncombustible environment. It is NOT designed for determining leaks in a combustible environment. A combustible gas detector should be used for determining combustible gas leaks under any condition. The Restek Electronic Leak Detector may be used for determining trace amounts of hydrogen in a GC environment only.

**Routine maintenance includes inspection of the probe tip, internal/external tubing and a battery replacement.

**i tech tip**

Avoid using liquid leak detectors on a GC! Liquids can be drawn into the system.

for more info

See **page 273** for more information on the Restek Electronic Leak Detector.