



Parker Balston® PEM Hydrogen Generators

- Proton Exchange Membrane (PEM) cell eliminates the need for liquid electrolytes.
- Reliably generate 99.9995% pure hydrogen—for better chromatography.
- Eliminates high-pressure cylinders—greater convenience and improved lab safety.
- Compact unit, requiring only one square foot of bench space.
- Quick and easy to service and maintain; unique display lighting changes color for easy status checks and water level indication.
- Comes with a set of universal power adapters for US, European, and Asian plug types.



Fuel-grade high purity hydrogen generators are safer alternatives to high-pressure gas cylinders. The new Proton Exchange Membrane (PEM) cell eliminates the use of liquid electrolytes with hydrogen generators. Deionized water is all that is required to generate hydrogen for weeks of continuous operation. With an output capacity of up to 510 cc/minute, one generator can supply 99.9995% pure hydrogen for up to several FIDs. Based on cylinder gas savings alone, a hydrogen generator pays for itself in one or two years.

Produced and supported by an ISO 9001 registered organization, Parker Balston® hydrogen generators are the first built to meet the toughest laboratory standards in the world: CSA, UL, cUL, and CE Mark. A great safety feature is the built-in sensing circuit, which shuts the generator down if a hydrogen leak is detected.

Specifications

Purity:	99.9995% pure hydrogen
Delivery Pressure:	10-100 psig \pm 1 psig (69-689 kPa \pm 7 kPa)
Outlet Port:	$\frac{1}{8}$ " compression
Electrical Requirements:	100-230 VAC/50-60 Hz
Physical Dimensions:	17.12" h x 13.46" w x 17.95" d (43.48 x 34.19 x 45.6 cm)
Shipping Weight:	40 lbs. (18 kg) dry



free literature

Parker Balston® Hydrogen Generators


Download your free copy from www.restek.com
lit. cat.# 580053A

Description	Model #	Capacity	qty.	cat.#	price
Hydrogen Generator	H2PEM-100	100cc/min.	ea.	23065	
Hydrogen Generator	H2PEM-165	165cc/min.	ea.	23066	
Hydrogen Generator	H2PEM-260	260cc/min.	ea.	23067	
Hydrogen Generator	H2PEM-510	510cc/min.	ea.	23068	

Replacement and Maintenance Components for Hydrogen Generators (for all models listed above)

Replacement Desiccant Cartridge for H2PEM Generators	ea.	23069
6-Month Maintenance Kit for H2PEM Generators Includes: 1 deionizer cartridge, 1 water filter, 3 environmental filters	kit	23070
24-Month Maintenance Kit for H2PEM Generators Includes: 1 deionizer cartridge, 1 water filter, 3 environmental filters, 1 water level sensor, 1 water pump, and 1 desiccant cartridge	kit	23071

The combination of high-purity gas and gas purifying traps can save analytical time in the long run.

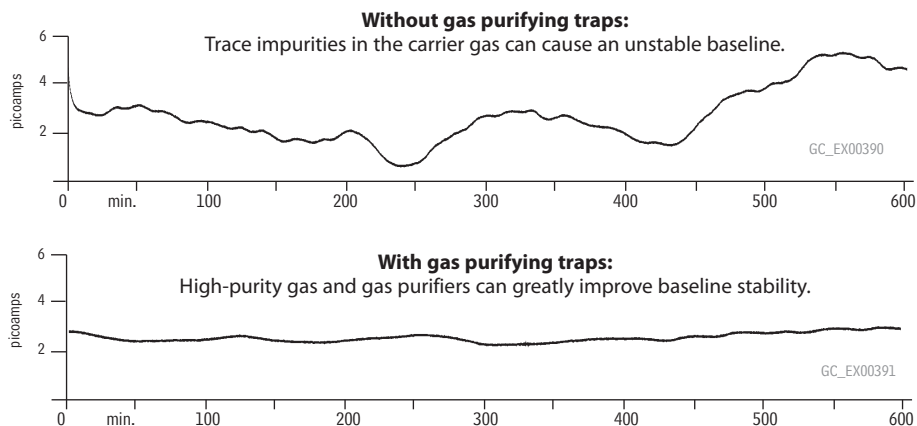


Using purifying gas traps in combination with a gas generator can improve chromatography tremendously and improve column lifetime.

Parker Balston® Model FID-1000 and FID-2500 Gas Stations

- Single unit produces UHP zero air from house compressed air and 99.9995% pure hydrogen from deionized water.
- Ideal for supplying up to 5-6 FIDs.
- Eliminates inconvenient and dangerous gas cylinders.
- Silent operation, minimal operator attention required.

Parker Balston® Gas Stations provide both UHP grade hydrogen gas and zero grade air for flame ionization detectors. The system is specifically designed to supply gas to FIDs and to support flame thermionic and flame photometric detectors. The units produce zero air by purifying compressed air to a total hydrocarbon concentration of 0.1 ppm or less (measured as methane).

The hydrogen generators produce hydrogen gas from deionized water, using the principle of electrolytic dissociation of water and hydrogen proton conduction through a proton exchange membrane cell.

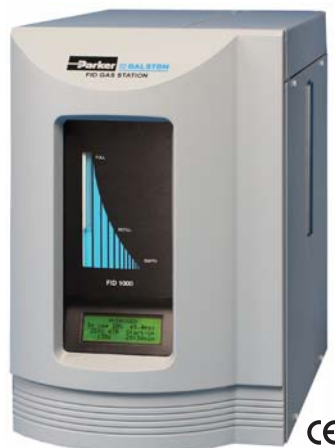

Produce zero air and pure hydrogen from one unit!

Specifications

Hydrogen Purity:	99.9995%	Power:	120 VAC/amp, 60 Hz, 400 watts
Zero Air Purity:	FID-1000:	Hydrogen Outlet Pressure:	60 psig (414 kPa)
	< 0.1 ppm total hydrocarbons as methane	Zero Air Outlet Pressure:	40-125 psig* (276-862 kPa)
	FID-2500:	Inlet Connection:	1/4" NPT (female)
	< 0.05 ppm total hydrocarbons as methane	Outlet:	1/8" compression
Max. Hydrogen Flow Rate:	FID-1000: 90 cc/min.	Dimensions:	16.5"h x 10.5"w x 17"d
	FID-2500: 250 cc/min.		(42 cm x 27 cm x 43 cm)
Max. Zero Air Flow Rate:	FID-1000: 1000 cc/min.	Weight:	53 lbs. (24kg)
	FID-2500: 2500 cc/min.		

*Zero air inlet requires minimum of 40 psig (276 kPa) compressed air pressure.



free literature

FID Gas Stations


Download your free copy from www.restek.com
lit. cat.# 580051

ordering note

For **international orders**, please add the appropriate power cord suffix from the table below.

Select the power cord you need.

Compare baselines produced by a Parker Balston® FID Gas Station and bottled fuel air. The baseline produced by the Parker Balston® Generator is flat, with no fluctuations or peaks; the chromatogram from the bottled air fuel supply has many peaks ranging from 0.25 ppm to -0.25 ppm total hydrocarbons.

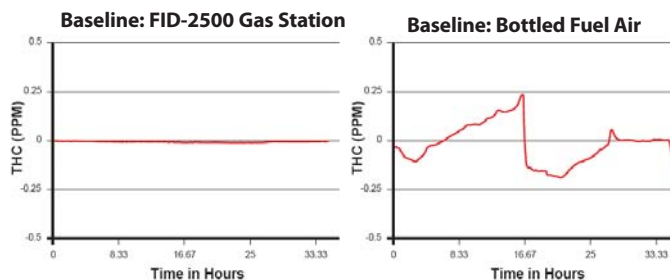

Description	Model #	qty.	cat. #	price
Gas Station	Model FID-1000 (ideal for 1-2 FIDs)	ea.	20177	
Gas Station	Model FID-2500 (ideal for 5-6 FIDs)	ea.	24913	

Replacement Components for FID Gas Stations

Resin Bed Cartridge for Hydrogen Generators in FID-1000 and FID-2500 Gas Stations		ea.	24914	
Replacement Desiccant Cartridge		ea.	21671	
FID Gas Station Maintenance Kit				
Includes: 1 desiccant cartridge, 1 resin bed cartridge, 1 filter cartridge		ea.	24915	

International Power Cord Sets

Just add the proper suffix to the catalog number for the gas generator you are ordering.

Location	qty.	cat.# suffix	price
United Kingdom (230VAC, 50/50Hz)	ea.	-550	
European (230VAC, 50/60Hz)	ea.	-551	
IEC Connector Only (230VAC, 50/60Hz)	ea.	-552	
Japanese (200VAC, 50/60Hz)	ea.	-556	
Japanese for Zero Air (100VAC, 50/60Hz)	ea.	-553	
Japanese for Hydrogen (100VAC, 50/60Hz)	ea.	-554	
Japanese for Nitrogen (100VAC, 50/60Hz)	ea.	-555	



www.restek.com 305



21654



Parker Balston® Nitrogen Gas Generators

- Produces ultra-pure nitrogen (up to 99.9995%).
- Flows from 1 to 75+ lpm.
- Require only a compressed air source and 110 volt AC power.
- Typical applications include GC carrier gas, make-up gas, and low flow sample concentrators.
- Maintenance kits include replacement filters.

Specifications

	Model HPN2-1100 or UHPN2-1100	Model HPN2-2000
Maximum Nitrogen Flow Rate:	See Flow Table	2 lpm
Nitrogen Purity:	99.9995%	99.99%
Minimum/Maximum Inlet Pressure:	60 psig/125 psig (414/862kPa)	75 psig/120 psig (517/827kPa)
Electrical Requirements:	120 VAC/60 Hz	120 VAC/60 Hz
Dimensions:	35" h x 12" w x 16" d (89cm x 30cm x 41cm)	35" h x 12" w x 16" d (89cm x 30cm x 41cm)
Shipping Weight:	115 lbs. (52 kg)	115 lbs. (52 kg)

Flow Table for Models HPN2-2000, HPN2-1100, and UHPN2-1100

Inlet Air Pressure	Maximum Outlet Flow (cc/min.)	
	Models HPN2-1100 and UHPN2-1100	Maximum Outlet Pressure
125 psig (862kPa)	1100	85 psig (586kPa)
110 psig (758kPa)	1000	75 psig (517kPa)
100 psig (689kPa)	900	65 psig (448kPa)
90 psig (621kPa)	800	60 psig (414kPa)
80 psig (552kPa)	700	50 psig (345kPa)
70 psig (483kPa)	600	45 psig (310kPa)
60 psig (414kPa)	500	35 psig (241kPa)
Model HPN2-2000		
75-120 psig (517-827kPa)	2000	90 psig (621kPa)



free literature

Nitrogen Generators

Download your free copy from www.restek.com

lit. cat.# 580052

Note:

Models HPN2-1100 and HPN2-2000 do not remove hydrocarbons.

*Power consumption is:

Model HPN2-1100

= 25 Watts

Model UHPN2-1100

= 700 Watts

Model HPN2-2000

= 25 Watts

ordering note

For **international orders**, please add the appropriate power cord suffix from the table below.

Nitrogen Generators	Model #	qty.	cat.#	price
Nitrogen Generator	HPN2-2000 (high purity)	ea.	21654	
Nitrogen Generator	HPN2-1100 (ultra-high purity)	ea.	21653	
Nitrogen Generator	HPN2-1100 with European Cord Set	ea.	21653-551	
Nitrogen Generator	HPN2-1100 with IEC Connector Only	ea.	21653-552	
Nitrogen Generator	UHPN2-1100 (ultra-high purity zero grade)	ea.	20697	
Maintenance Kits	Model #	qty.	cat.#	price
Maintenance Kit	for HPN2-1100, HPN2-2000, 76-96, 76-92	kit	21649	
Maintenance Kit	for UHPN2-1100, 76-94	kit	21655	

International Power Cord Sets

Just add the proper suffix to the catalog number for the gas generator you are ordering.

Location	qty.	cat.# suffix	price
United Kingdom (230VAC, 50/50Hz)	ea.	-550	
European (230VAC, 50/60Hz)	ea.	-551	
IEC Connector Only (230VAC, 50/60Hz)	ea.	-552	
Japanese (200VAC, 50/60Hz)	ea.	-556	
Japanese for Zero Air (100VAC, 50/60Hz)	ea.	-553	
Japanese for Hydrogen (100VAC, 50/60Hz)	ea.	-554	
Japanese for Nitrogen (100VAC, 50/60Hz)	ea.	-555	

Select the power cord you need.

Parker Balston® Nitrogen Gas Generators for LC/MS

- Turn compressed air into ultra-pure nitrogen (up to 99.5%).
- Flows from 1 to 44 lpm.
- Models N2-04, N2-14, N2-22, and N2-35 require no electricity.
- Safe, reliable, low maintenance.
- Maintenance kits include replacement filters.



22130



22129



Specifications	NitroFlow Lab	N2-04	N2-14 or N2-14A	N2-22 or N2-22A	N2-35 or N2-35A
Maximum Nitrogen Flow Rate:	32 lpm	8 lpm	78scfh at 95% purity	N2-22: 44 lpm N2-22A: 29 lpm	44 lpm
Nitrogen Purity:	99.50%	99%	95.0%–99.5%	99%	99%
Min/Max Inlet Pressure:	N/A	60 psig/145 psig	60 psig/145 psig	60 psig/145 psig	60 psig/145 psig
Electrical Requirements:	120VAC/60Hz	None	N2-14: None N2-14A: 120VAC/60Hz	N2-22: None N2-22A: 120VAC/60Hz	N2-35: None N2-35A: 120VAC/60Hz
Dimensions:	27.6" h x 35.4" w x 12.2" d (70cm x 90cm x 31cm)	11" h x 13" w x 16" d (27cm x 34cm x 41cm)	50" h x 16" w x 16" d (127cm x 41cm x 41cm)	50" h x 16" w x 16" d (127cm x 41cm x 41cm)	50" h x 16" w x 16" d (127cm x 41cm x 41cm)
Shipping Weight:	205 lbs. (93 kg)	43 lbs. (20 kg)	N2-14: 75 lbs. (34 kg) N2-14A: 80 lbs. (36 kg)	N2-22: 101 lbs. (46 kg) N2-22A: 106 lbs. (48 kg)	N2-35: 115 lbs. (52 kg) N2-35A: 119 lbs. (54 kg)

Nitrogen Generators for LC/MS Model #

	Model #	qty.	cat.#	price
Nitrogen Generator for LC/MS	NitroFlow Lab Model, 32 lpm max. flow	ea.	22129	
Nitrogen Generator for LC/MS	N2-04 Model for ELSD, 8 lpm max. flow	ea.	22130	
Nitrogen Generator for LC/MS	N2-14 (general purpose) 78 scfh max. flow at 95% purity	ea.	20677	
Nitrogen Generator for LC/MS	N2-14 with European Power Cord	ea.	20677-551	
Nitrogen Generator for LC/MS	N2-14A (general purpose w/oxygen analyzer) 78 scfh max. flow at 95% purity	ea.	21652	
Nitrogen Generator for LC/MS	N2-22 Model, 44 lpm max. flow	ea.	22131	
Nitrogen Generator for LC/MS	N2-22A Model, 29 lpm max. flow	ea.	22132	
Nitrogen Generator for LC/MS	N2-35 Model, 44 lpm max. flow	ea.	22133	
Nitrogen Generator for LC/MS	N2-35A Model, 44 lpm max. flow	ea.	22134	

Maintenance Kits

	Model #	qty.	cat.#	price
Maintenance Kit	for N2-14, N2-14A, 75-72, 75-720NA	kit	21648	
Maintenance Kit with Carbon Filter	for N2-14, N2-14A, 75-72, 75-720NA	ea.	22135	

ordering note

For **international orders**, please add the appropriate power cord suffix from the table on the previous page.

Restek Electronic Leak Detector...

Go to www.restek.com/leakdetector for details.

...and Introducing the NEW Restek ProFLOW 6000 Electronic Flowmeter

Go to www.restek.com/flowmeter for details.

RESTEK

CHROMALYTIC +61(0)3 9762 2034
ECHnology Pty Ltd

Australian Distributors
Importers & Manufacturers
www.chromtech.net.au

11/12

www.restek.com **307**



Parker Balston® Zero Air Generators

- Turn in-house compressed air into ultra-pure air (<0.1 ppm total hydrocarbons).
- Remove hydrocarbons to less than 0.1 ppm by catalytic oxidation.
- Operate at 40 to 125 psi (276-862 kPa).
- Typical payback is less than one year, based on cylinder costs.
- Install easily and take up little bench space.
- Maintenance kits include a one year supply of prefilters and final filter.



Specifications

Maximum Zero Air Flow Rate:	75-83NA	1 lpm
	HPZA-3500	3.5 lpm
	HPZA-7000	7 lpm
	HPZA-18000	18 lpm
	HPZA-30000	30 lpm
Outlet Hydrocarbon Concentration (as methane):	75-83NA	< 0.1 ppm
	HPZA-30000	< 0.1 ppm
	Other Models	< .05 ppm
Minimum/Maximum Inlet Air Pressure:	40 psig/125 psig (276/862 kPa)	
Maximum Inlet Hydrocarbon Concentration (as methane):	100 ppm	
Pressure Drop at Maximum Flow Rate:	4 psi (28 kPa) differential	
Maximum Inlet Air Temperature:	78°F (25°C)	
Inlet/Outlet Ports:	1/4" NPT (female)	
Start-up Time to Specified Hydrocarbon Concentration:	45 minutes	
Electrical Requirements:	75-83NA	120 VAC/60 Hz, 0.5 amps
	Other Models	120 VAC/60 Hz, 3.5 amps
Dimensions:	75-83NA	12"h x 10"w x 3"d (30 cm x 25 cm x 8 cm)
	Other Models	16"h x 11"w x 13"d (42 cm x 27 cm x 3 cm)
Shipping Weight:	75-83NA	7 lbs. (3 kg)
	Other Models	41 lbs. (19 kg)

Model	Number of FIDs*
75-83NA	Up to 3
HPZA-3500	Up to 11
HPZA-7000	Up to 23
HPZA-18000	Up to 60
HPZA-30000	Up to 100

*based on a 300 cc/min. fuel air rate

Zero Air Generator	Model #	Capacity	qty.	cat. #	price
Zero Air Generator	75-83NA	1000cc/min.	ea.	20684	
Zero Air Generator	75-83NA with United Kingdom Power Cord	1000cc/min.	ea.	20684-550	
Zero Air Generator	HPZA-3500	3500cc/min.	ea.	20680	
Zero Air Generator	HPZA-3500 with European Power Cord	3500cc/min.	ea.	20680-551	
Zero Air Generator	HPZA-7000	7000cc/min.	ea.	20681	
Zero Air Generator	HPZA-18000	18,000cc/min.	ea.	20682	
Zero Air Generator	HPZA-30000	30,000cc/min.	ea.	20683	

Maintenance Kits (includes a one-year supply of prefilters and final filter)	Model #	Capacity	qty.	cat. #	price
Maintenance Kit	for 75-83NA		kit	21646	
Maintenance Kit	for HPZA-3500, HPZA-7000, HPZA-18000, HPZA-30000		kit	21647	

Replacement Catalyst Towers	Model #	Capacity	qty.	cat. #	price
Replacement Catalyst Tower	for 75-83NA	1000cc/min.	ea.	22005	
Replacement Catalyst Tower	for HPZA-3500	3500cc/min.	ea.	22004	
Replacement Catalyst Tower	for HPZA-7000	7000cc/min.	ea.	22006	
Replacement Catalyst Tower	for HPZA-18000	18,000cc/min.	ea.	22007	
Replacement Catalyst Tower	for HPZA-30000	30,000cc/min.	ea.	22008	

ordering note

For **international orders**, please add the appropriate power cord suffix from the table below.



free literature

Zero Air Generators

Download your free copy from www.restek.com
lit. cat.# 580050

International Power Cord Sets

Just add the proper suffix to the catalog number for the gas generator you are ordering.

Location	qty.	cat.# suffix	price
United Kingdom (230VAC, 50/50Hz)	ea.	-550	
European (230VAC, 50/60Hz)	ea.	-551	
IEC Connector Only (230VAC, 50/60Hz)	ea.	-552	
Japanese (200VAC, 50/60Hz)	ea.	-556	
Japanese for Zero Air (100VAC, 50/60Hz)	ea.	-553	
Japanese for Hydrogen (100VAC, 50/60Hz)	ea.	-554	
Japanese for Nitrogen (100VAC, 50/60Hz)	ea.	-555	