


BHA and BHT
Rtx®-50

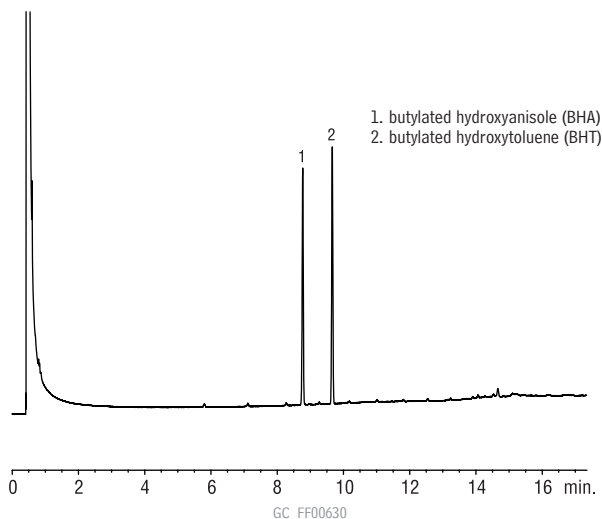

Column: Rtx®-50, 30m, 0.53mm ID, 0.50 μ m (cat.# 10540)
Sample: 50ppm each in methanol
Inj.: 1.0 μ L direct, gooseneck splitless inlet liner, 4mm (cat.# 20798)
Inj. temp.: 280°C
Carrier gas: helium, constant pressure
Linear velocity: 60cm/sec. @ 50°C
Oven temp.: 50°C to 240°C @15°C/min. (hold 3 min.)
Det.: FID @ 280°C

Chromatogram Search Tool

Search by compound name, synonym, CAS # or keyword

www.restek.com/chromatograms

Tocopherols and Tocotrienols
Rtx®-20

