


Vitamin A Compounds  
Ultra II® Aromax



Peak List	Ret. Time (min.)
1. retinol	0.832
2. β-carotene	0.950

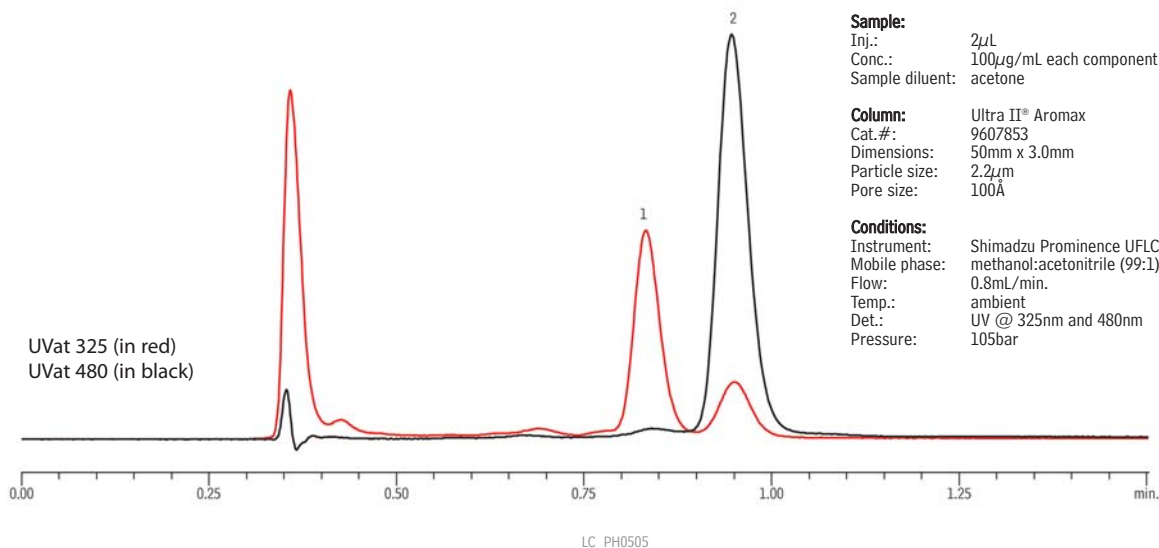

**Sample:**  
Inj.: 2µL  
Conc.: 100µg/mL each component  
Sample diluent: acetone

**Column:** Ultra II® Aromax  
Cat.#: 9607853  
Dimensions: 50mm x 3.0mm  
Particle size: 2.2µm  
Pore size: 100Å

**Conditions:**  
Instrument: Shimadzu Prominence UFLCXR  
Mobile phase: methanol:acetonitrile (99:1)  
Flow: 0.8mL/min.  
Temp.: ambient  
Det.: UV @ 325nm and 480nm  
Pressure: 105bar

## Chromatogram Search Tool

Search by compound name, synonym,  
CAS # or keyword

[www.restek.com/chromatograms](http://www.restek.com/chromatograms)

