


Gas Line Installation Kit

The Gas Line Installation Kit includes everything you will need to connect a single gas cylinder to an SRI GC.

- **Cylinder Pressure Regulator with 0-100psi Output**
- **50-foot length of 1/8" Copper tubing**
- **Stainless steel Gas Line Filter**
- **Extra Swagelok nuts & ferrules**
- **Handy tubing cutter**

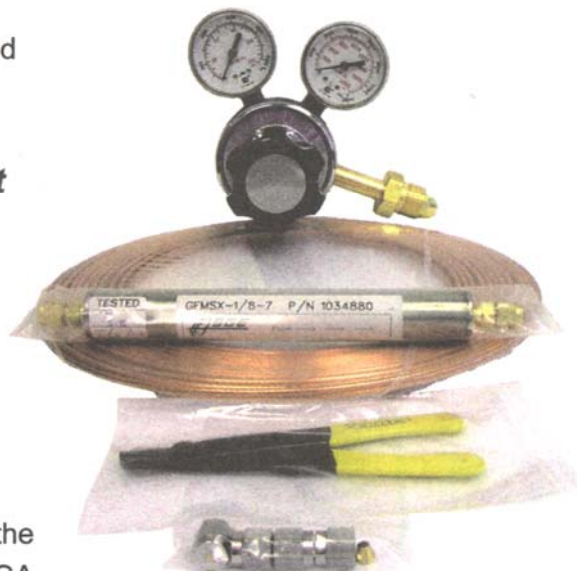

Each type of gas cylinder has a different type of CGA connection to the cylinder. Air is typically CGA 590 or 346, helium and nitrogen are CGA 580, while hydrogen and argon-methane (P5) are CGA 350. Specify the CGA type when ordering. Each cylinder regulator is supplied with a 1/8" Swagelok outlet fitting for easy connection to the copper tubing, and the hydrogen regulator is also equipped with a flow restrictor to limit the escape of gas in the event of a leak. The air kit is supplied with both CGA 590 and CGA 346 inlet fittings.

8600-C590	Air Gas Line Kit	\$ 395.00
8600-C580	Helium/Nitrogen Gas Line Kit	\$ 395.00
8600-C350	Hydrogen/Argon-Methane Gas Line Kit	\$ 395.00

GC Maintenance Kit

The GC Maintenance Kit includes most parts which could fail, and is designed especially for our export customers who may have difficulty returning the GC to the factory for service.

- **Replacement chips for most electronic circuits**
- **Replacement heaters for column oven & detectors**
- **Assorted graphite & Vespel ferrules**
- **Assorted Swagelok nuts & ferrules**
- **Digital voltmeter for troubleshooting**
- **Wide-bore column adapter**
- **Type K thermocouple**
- **EPC pressure sensor**

8600-MAIN	GC Maintenance Kit (specify 120 or 220 VAC)
------------------	--