


PeakSimple Chromatography Data Systems



4 channel data system
built into an 8610C GC

- **Built into Every SRI HPLC, 8610 and 310 GC**
- **Available separately for use with most GCs and HPLCs**
- **Easy Connection to your Windows™ PC**
- **1, 4, and 6 Channel Models available**
- **PeakSimple Software Included**
- **Serial port or USB models**

4 channel stand-alone
data system

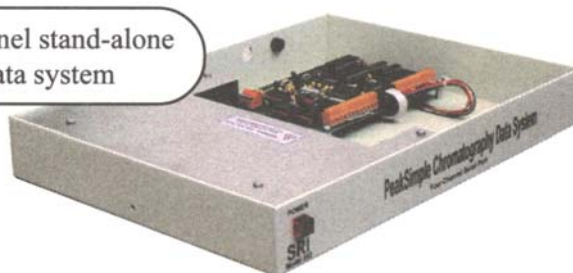

PeakSimple Chromatography Data Systems consist of hardware and software. The hardware is available as a stand-alone data system for connection to almost any model GC, HPLC or CE system. The same hardware is supplied as standard equipment with every SRI HPLC, 8610 and 310 GC. No hardware is installed in your computer, so a portable laptop may be used instead of a full-sized desktop PC. PeakSimple chromatography acquisition and integration software for Windows is provided with each data system, and updates are FREE from the SRI website: www.srigc.com. The data system hardware comes in three versions:

- 1) Single channel Model 203 for one detector
(RS232 serial connection to PC or USB with optional adaptor)
- 2) 4 channel serial Model 202 for up to 4 detectors on 1 or 2 instruments
(RS232 serial connection to PC or USB with optional adaptor)
- 3) 6 channel USB Model 302 for up to 6 detectors on 1-4 instruments
(USB connection to PC)

All three models use the same PeakSimple software.

