


Model: STE-8L(N)

# AUTOCLAVE

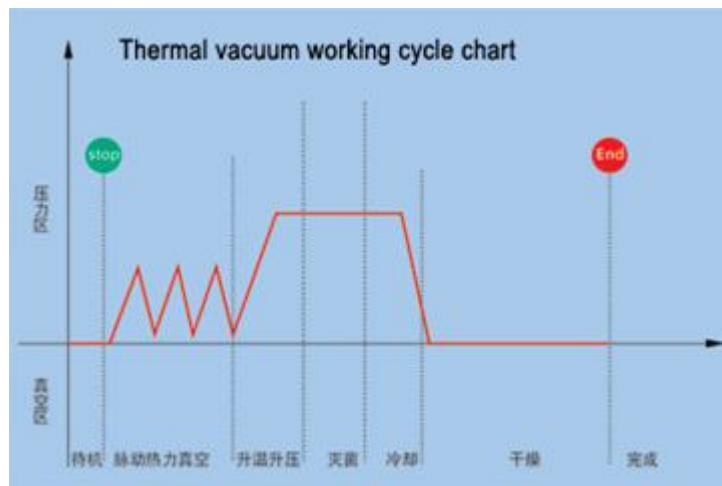

**Description:**

Adopting European N standard, with thermal vacuuming and drying function, it can be used for A type hollow instruments and solid instruments

**Class N sterilization working cycle chart:**

**Technical feature:**

- \* Sterilization temperature: 121, 134 Cent degrees
- \* Drying Program: Thermal Vacuum drying
- \* Digital display
- \* Sterilization record: mini external printer is available
- \* Safety: safety valve, temperature control, pressure protection system, temperature protector of chamber wall.
- \* Sterilization and drying time adjustable help you to get the best sterilization effect.
- \* Open type water storage tank
- \* Quick drain connector

**Specifcation:**

Type	Voltage (V/Hz)	Power (W)	Capacity (L)	Chamber (mm)	Outside size (mm)	Package Size (mm)	Net Weight (KG)	Gross Weight (KG)
STE-8L(N)	AC220V/50Hz	1100W	8L	Φ170×320	345×530×330	640×460×440	29Kg	34Kg

Temperature (°C)	Pressure (bar)	Holding time (min)	Total time (min)	TYPE	MAXLOAD (kg)
134	2.10	6	30~40	Unwrapped hollow material	5.00
121	1.10	30	35~50	Single-wrapped solid material	4.00