

PEEK no-twist “one-piece” fittings

10-32

- Snap-in ferrule rotates freely
- Choice of ferrule materials
- Choice of fitting lengths

No-twist fittings offer the convenience of a one-piece fitting while solving a problem inherent to such designs. In other one-piece designs, the ferrule rotates against the fitting detail, creating particulates. The no twist design has a separate ferrule that snaps into the nut, so it's attached but still free to avoid rotation during tightening.



Since the ferrule is not machined onto the nut, it can be made from a different material. We offer several PEEK nuts, with ferrules of PEEK, CTFE, or glass-filled PEEK.

Sold in packages of 5.


Ferrule material:	PEEK	Glass-filled PEEK	CTFE
Description	Product No.	Product No.	Product No.
Fingertight	JR-ZNF1FPK-5	JR-ZNF1FPKG-5	JR-ZNF1FKF-5
Hex, short	JR-ZNF1PK-5	JR-ZNF1PKG-5	JR-ZNF1KF-5
Hex, medium	JR-MZNF1PK-5	JR-MZNF1PKG-5	JR-MZNF1KF-5
Hex, long	JR-LZNF1PK-5	JR-LZNF1PKG-5	JR-LZNF1KF-5
Replacement ferrules	JR-ZNF1PK-F-5	JR-ZNF1PKG-F-5	JR-ZNF1KF-F-5

SPECS


- For 1/16" OD tubing
- Materials:
Nut: PEEK
Ferrule: PEEK, glass-filled PEEK, or CTFE
- 10-32 threads
- Pressure rating:
<350 bar (5000 psi)
(Varies with tubing material and ID)

OPTIONS


- Optional ferrule materials available: FEP, PFA, PTFE, and glass-filled PTFE. Contact us for availability.
- Some 1/32" versions are available. Call for details.

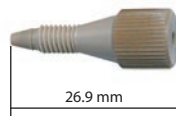


Fingertight

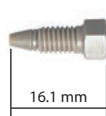

Hex, short

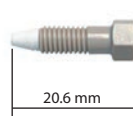

Hex, medium

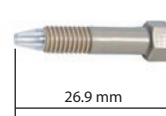

Hex, long

RELATED PRODUCTS

- 1/16" OD tubing for high pressure applications 5-8
- PEEK starter kit 90
- ValvTool for tightening hex-head fittings 88
- Color-It fingertight adapters. 17

