

PEEK ferrules

NEW

UHPLC

SPECS

- For 1/16" OD tubing
- Material: PEEK, special black
- Pressure rating on SS tubing and with SS nut: (preliminary) 1170 bar (17,000 psi)

- Biocompatible 1/16" ferrule for ultra-high pressures
- Alternative to titanium ferrules in applications with titanium valves
- Excellent in combination with stainless nuts

Sold in packages of 5.

Product No.
JR-5010-5



PEEK nuts with single ferrules

10-32

SPECS

- Material: PEEK
- Pressure rating: <280 bar (4000 psi) (Varies with tubing material and ID)


- Reuseable nut, replaceable ferrule
- For 1/32", 1/16", and 1/8" tubing
- Suitable for pressures to 280 bar (4000 psi)

Fingertight nuts have a knurled surface designed to allow sufficient sealing force to be applied to the ferrule without wrenches. Hex style nuts allow wrench tightening; however, since they are polymeric, they can break, and are recommended for use only when space is limited and fingers will not fit.

Sold in packages of 5.



	Tubing OD	Description	Thread	Product No.
Nuts	1/32"	Fingertight	6-40	JR-ZN.5FPK-5
	1/16"	Fingertight	10-32	JR-5570-5
	1/8"	Hex-head	5/16-24	JR-ZN2PK-5
Single ferrules	1/32"			JR-ZF.5PK-5
	1/16"			JR-5003-5
	1/8"			JR-ZF2PK-5

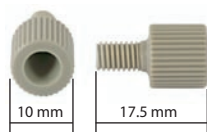

1/16" fingertight



1/8" hex-head

MORE INFORMATION

About PEEK 99
Chemical resistance 96-98

RELATED PRODUCTS

1/16" OD tubing for high pressure applications 5-8
PEEK starter kit 90
PEEK tubing 5-8
Stainless steel nuts ... 31
ValvTool for tightening hex-head fittings 88
Color-It fingertight adapters 17

