

PEEK fingertight nuts (flangeless)

1/4-28

SPECS

- For 1/16" and 1/8" OD tubing
- Materials
Nut: PPS
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)

- Compatible with all 1/4-28 flat-bottom ports
- Excellent chemical resistance
- For 1/16" or 1/8" OD tubing
- True fingertight – no tools required




PEEK flangeless nuts in combination with our inverted ferrule do not require a flanging tool. A hand-tightening procedure makes leak-free seals.

With this fitting system, the inert ETFE ferrule is the only part besides the tubing which comes into contact with the fluid stream – the nut is never in contact.

Ferrules are not included. Order separately.



Tubing OD:		1/16"		1/8"	
Description	Material	Qty/pkg	Product No.	Product No.	
Nut for use with inverted ferrule	PEEK	10	JR-20124-10	JR-20116-10	
Inverted ferrules	ETFE	10	JR-041-10	JR-051-10	

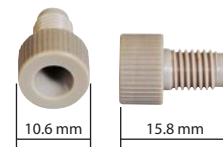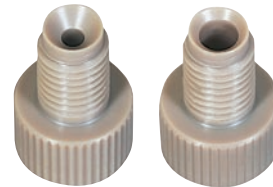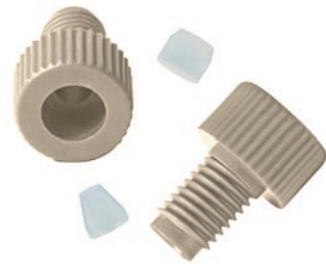

Color-coded fingertight nuts (flangeless)

1/4-28

SPECS

- For 1/16" and 1/8" OD tubing
- Materials
Nut: Polyacetal
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:

OD	Bar	Psi
1/16"	70	1000
1/8"	35	500


 (Varies with tubing material and ID)

- Polyacetal – economical alternative to PEEK
- Compatible with all 1/4-28 flat-bottom ports
- For 1/16" or 1/8" OD tubing
- True fingertight – no tools required

Polyacetal flangeless nuts in combination with the inverted ferrule do not require a flanging tool. A hand-tightening procedure makes leak-free seals. With this fitting system, the inert ETFE ferrule is the only part besides the tubing which comes into contact with the fluid stream – the nut is never in contact. Six colors are available.



Tubing OD:		1/16"		1/8"	
Description	Color	Qty/pkg	Product No.	Product No.	
Nuts	Green	10	JR-040GR-10	JR-050GR-10	
	Blue	10	JR-040BL-10	JR-050BL-10	
	Yellow	10	JR-040YL-10	JR-050YL-10	
	Red	10	JR-040RE-10	JR-050RE-10	
	White	10	JR-040WH-10	JR-050WH-10	
	Black	10	JR-040BK-10	JR-050BK-10	
1" long nuts	White	10	JR-040WL-10	JR-050WL-10	
Inverted ETFE ferrules		10	JR-041-10	JR-051-10	

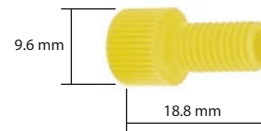

TECH TIP

The inert ETFE ferrule used with these nuts is the only part besides the tubing which is in contact with the fluid stream.

