




Trapping/flow-through flowpath – STF configuration

The STF selector is a variation of the ST flowpath, with the single difference that the non-selected streams are returned to their own vents or sources rather than being dead-ended or trapped as they are in the standard ST configuration. For an application suggestion, see page 138.

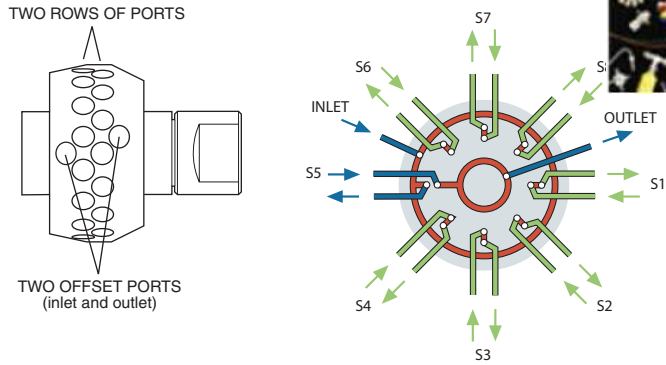

1/16" fittings, 0.75 mm ports (.030")

MWType

Low pressure	Includes 2" standoff. Ask about closemount assembly if valve will not be heated.	Standard electric actuators: 110 VAC for USA 110/230 VAC to 24 VDC power supply for international Microelectric actuators: 24 VDC (includes a 110/230 VAC to 24 VDC power supply)	SPECS 200 psi gas 200°C max Nitronic 60 body Valcon E rotor	
STF Trap/ flow-throw				
1/16"				0.75 mm

OPTIONS

- 4 and 8 positions available
- 3", 4", and 6" standoffs
- Materials: Hastelloy C, Inconel 600, Monel 400, Nickel 200, Nitronic 50, Titanium (see pages 254-255)

	6 Position		10 Position		12 Position		16 Position	
	Prod No	Price	Prod No	Price	Prod No	Price	Prod No	Price
Manual (not recommended)	2CSTF6MWE		2CSTF10MWE		2CSTF12MWE		2CSTF16MWE	
With air actuator	A2CSTF6MWE		A2CSTF10MWE		A2CSTF12MWE		A2CSTF16MWE	
With standard elec actuator	E2CSTF6MWE		E2CSTF10MWE		E2CSTF12MWE		E2CSTF16MWE	
With microelectric actuator	EMT2CSTF6MWE		EMT2CSTF10MWE		EMT2CSTF12MWE		EMT2CSTF16MWE	
Replacement valve	DCSTF6MWE		DCSTF10MWE		DCSTF12MWE		DCSTF16MWE	
Replacement rotor	SSACSTF6MWE		SSACSTF10MWE		SSACSTF12MWE		SSACSTF16MWE	

MORE INFORMATION

Application..... page 138

Actuators

Air 194

Microelectric .. 190-191

Standard elec..... 193

Materials

Metals..... 254-255

Polymers 256

Valve rotors..... 257

Mounting hardware

Closemount 208

Standoff..... 205



1/8" fittings, 1.0 mm ports (.040")

MW Type

SPECS
 200 psi gas
 200°C max
 Nitronic 60 body
 Valcon E rotor

Includes 2" standoff.
 Ask about closemount
 assembly if valve will
 not be heated.

Standard electric actuators:
 110 VAC for USA
 110/230 VAC to 24 VDC power supply for international
 Microelectric actuators:
 24 VDC (includes a 110/230 VAC to 24 VDC power supply)

Low pressure
 STF
 Trap/ flow-throw
 1/8" 1.0 mm

OPTIONS

- 4 and 8 positions available
- 3", 4", and 6" standoffs
- Materials: Hastelloy C, Inconel 600, Monel 400, Nickel 200, Nitronic 50, Titanium (see pages 254-255)
- Larger bore available except 16 position

	6 Position		10 Position		12 Position		16 Position	
	Prod No	Price	Prod No	Price	Prod No	Price	Prod No	Price
Manual (not recommended)	2STF6MWE		2STF10MWE		2STF12MWE		2STF16MWE	
With air actuator	A2STF6MWE		A2STF10MWE		A2STF12MWE		A2STF16MWE	
With standard elec actuator	E2STF6MWE		E2STF10MWE		E2STF12MWE		E2STF16MWE	
With microelectric actuator	EMT2STF6MWE		EMT2STF10MWE		EMT2STF12MWE		EMT2STF16MWE	
Replacement valve	DSTF6MWE		DSTF10MWE		DSTF12MWE		DSTF16MWE	
Replacement rotor	SSASTF6MWE		SSASTF10MWE		SSASTF12MWE		SSASTF16MWE	

