

Stream Selectors



UHPLC and High Pressure Selectors

NEW Model C75 selectors offer pressure ratings of 15,000 psi and 10,000 psi with 1/32" and 1/16" fittings.

The **Model C5**, with Valco ZDV fitting details, is available with 4, 6, 8, or 10 positions. Stators are available in Nitronic 60 stainless, titanium, and Hastelloy C-22, with rotors of Valcon H, all of which are compatible with common HPLC solvents. PAEK stators are used in combination with Valcon E rotors.

The C5 valve is the backbone of the Cheminert **HPLC column selector system**, which includes two stream selection valves mounted on a single microelectric actuator. Columns are not included.

Model C75
pages 170 and 171


15,000 psi

10,000 psi



Model C5
page 172

5,000 psi

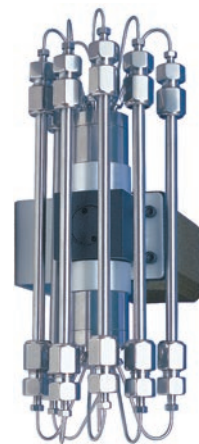

HPLC column selector system
page 173

5,000 psi

Column selector system

Low Pressure Selectors

With Valco Zero Dead Volume Fittings

Model C25Z valves have Valco ZDV fitting details, and are available in 4, 6, 8, 10, 12, and 14 position models.

Model C35Z valves have 1/16" Valco ZDV details, and are available in 20, 24, and 26 position models. This is a tapered rotor valve limited to 100 psi liquid. Rotors are made from Valcon E2, with valve body made from PPS.

With Cheminert 1/4-28 Fittings

The **Model C25** has female 1/4-28 threaded fitting details for direct connection of lines – no couplings are required. The C25 is available in 4, 6, 8, and 10 position models. Multicolored Cheminert 1/16" or 1/8" flangeless fittings are included. Order other fittings separately as required. Rotors are made of Valcon E2, a proprietary reinforced PTFE composite, with stators of PPS.

With Cheminert 1/2-20 Fittings

Model C45 valves feature 1/2-20 threaded fitting details for use with 1/4" OD tubing. This is a tapered rotor valve with large bore for high flow applications. Rotors are made from Valcon E2, with valve body made from PPS. Available in 4 and 6 port configurations.



Model C25Z
page 174

Low pressure


10-32 ZDV



Model C25
page 175

Low pressure


1/4-28 Internal

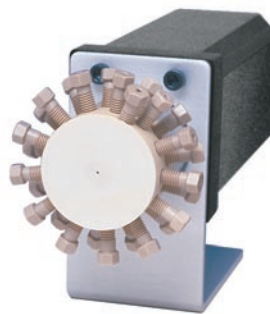

Model C35Z
page 176

Low pressure

10-32 ZDV

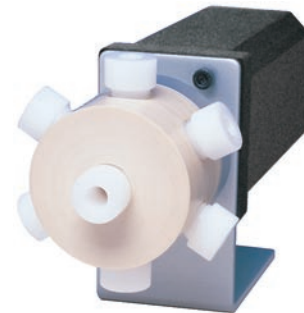

Model C45
page 177

Low pressure

1/2-20 Int.

CHROMalytic TECHNOlogy Pty Ltd AUSTRALIAN Distributors e-mail: sales@chromtech.net.au Tel: 03 9762 2034