



PEEK starter kit

In LC applications involving proteins, peptides, nucleic acids, or other samples of biological origin, metal systems may interact with samples or release transition metals that will deactivate columns. The PEEK starter kit facilitates replacement of stainless steel tubing, fittings, ferrules, mobile phase filters, etc., to create a biocompatible environment for samples and mobile phase.

Prod No

PEEK starter kit


JR-35P

Includes:

- 1 Plastic box
- 10 PEEK one-piece fittings, 10-32
- 5 PEEK handtight fittings
- 5 PEEK nuts, hex-head long
- 20 PEEK ferrules, double-ended 1/16"
- 1 PEEK union, HP body only, 10-32
- 2 Tubing elbows 90°
- 2 Tubing elbows 180°
- 1 PEEK filter, in-line,
incl. PAT frit 5 µm
- 1 Clean-cut tubing cutter
- 1 Last Drop PTFE filter 5 µm
- 3m PEEK tubing, 1/16" x 0.25 mm ID,
blue stripe
- 3m PEEK tubing, 1/16" x 0.50 mm ID,
orange stripe
- 1 Tweezers

Tightening tools for PEEK fittings

These handy tools make it fast and easy to tighten PEEK hex-head fittings. The red version is for use with the C360 series fittings shown on page 57. The green tool is for any 1/32" PEEK fitting with a 3/16" hex head nut, and the blue version fits the 1/4" hex common in fittings for 1/16" tubing.

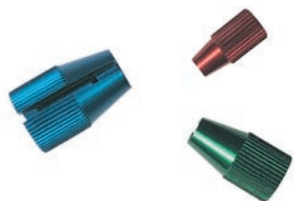

Color	For use with	Prod No
Red	360 µm fittings	C360ET
Green	1/32" fittings	CNFT
Blue	1/16" fittings	ZNFT

MORE INFORMATION

Hex-head PEEK fittings
 360 µm..... page 57
 1/32" 63-65
 1/16"
 High pressure... 63-65
 Low pressure 71

CHROMalytic TECHNOlogy Pty Ltd AUSTRALIAN Distributors e-mail: sales@chromtech.net.au Tel: 03 9762 2034