



No-Met biocompatible mobile phase filter

Stainless steel in the flowpath is not acceptable in a growing number of applications involving the separation of biomolecules. High salt buffers can corrode stainless steel, and the metal ions released from metallic filters may contaminate or otherwise react with the biomolecules of interest.

The No-Met polyethylene filter is designed for these applications, with inert polymeric fittings and 20 µm filter effectively eliminating metal contamination from the fluid path. Use them for IC and biochromatography applications.

Because they are hydrophobic, No-Met filters may initially require some priming with methanol or acetonitrile.



	<i>Prod No</i>	<i>Price</i>
No-met mobile phase filter, 1/8"	JR-32178	
Replacement element	JR-32179	

Stainless steel mobile phase filters and helium spargers

Mobile phase filters protect your HPLC system from small particles in the mobile phase. These filters are made from 316 stainless and PEEK or PTFE, and are suitable for use with most solvents.

Helium spargers offer an inexpensive way to prepare and maintain mobile phases free of dissolved gases. Connect these spargers to a regulated supply of helium gas (0-400 ml/min) to remove dissolved gases from the mobile phase. Spargers are made from 10 micron porosity stainless steel.




<i>Tubing OD</i>	<i>Porosity</i>	<i>Suggested Max. Flow Rate (ml/min)</i>	<i>Prod No</i>	<i>Price</i>
1/16"	2 µm	8	JR-367016-2	
1/16"	10 µm	20	JR-367016-10	
1/16"	20 µm	20	JR-367016-20	
1/8"	2 µm	8	JR-367008-2	
1/8"	10 µm	20	JR-367008-10	
1/8"	20 µm	20	JR-367008-20	

Mobile phase filters

Direct connect

Cheminert mobile phase filters provide point-of-use filtering of common HPLC or FIA solvents. They are designed to connect directly to 1/8" OD PTFE or PEEK tubing using a simple press fit. The filter housing is PTFE and includes a 2 or 10 micron titanium frit.

<i>Pore size</i>	<i>Prod No</i>	<i>Price</i>
2 micron	C-MPFT12	
10 micron	C-MPFT110	

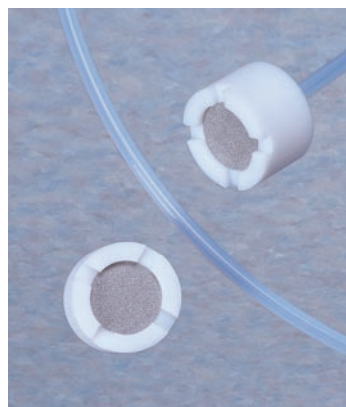

0.25 mm = .010"
0.50 mm = .020"
0.75 mm = .030"
1.0 mm = .040"
1.5 mm = .060"
2.0 mm = .080"
4.6 mm = .180"
6.0 mm = .236"
6.4 mm = .253"
7.0 mm = .275"
10.0 mm = .400"
27.0 mm = 1.08"
1/32" = 0.8 mm
1/16" = 1.6 mm
1/8" = 3.2 mm
1/4" = 6.4 mm
3/8" = 9.5 mm
1/2" = 12.7 mm