



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

- Tubing for chromatography
- High pressure fittings
 - Nanovolume
 - Stainless steel
 - Polymeric
- Low pressure fittings
- Filters & mobile phase filters
- Mobile phase handling systems
- Frits
- PEEK columns & guard columns
- Sample loops
- Accessories/miscellaneous

CUSTOMER SUPPORT

PEEK, polyacetal, and polypropylene color-coded nuts, fingertight adapters, unions, tees, crosses, and much more.

Low Pressure Polymeric Fittings

PPS nuts for flat bottom 1/4"-28 fitting details

Flangeless type, collapsible ferrule

Flangeless type, inverted ferrule

Flanged type

Color-Sleeve fingertight adapters

Color-coded fingertight nuts

PEEK, flangeless

Polyacetal, flangeless

Polypropylene, flanged

Flanged tubing connector

Unions

Bulkhead unions

Tees, crosses, and manifolds

Biocompatible mixing tee

Polypropylene adapters

Manual flow selection valves

Peristaltic tubing adapter set

Flanged fitting kit



TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous


CUSTOMER SUPPORT

Unions

- Easy connection of low pressure solvent lines
- Low dead volume
- Completely inert and biocompatible

These unions allow the connection of two pieces of 1/16" tubing, two pieces of 1/8" tubing, or one 1/8" to one 1/4" piece of tubing, quickly and reliably. The unions are inert and biocompatible.

VICI unions have 1/4"-28 flat bottom ports, so they can be used with collapsible ferrules, inverted ferrules, or flanged tube ends. Complete versions of the 0.50 mm and 0.75 mm bore unions are shipped with PPS flangeless nuts and ETFE ferrules for 1/16" tubing; 1.30 mm versions get PEEK nuts and ETFE ferrules for 1/8" tubing. The reducing union is supplied with polyacetal flangeless nuts, with an ETFE ferrule for the 1/8" connection and a CTFE ferrule for the 1/4" connection.

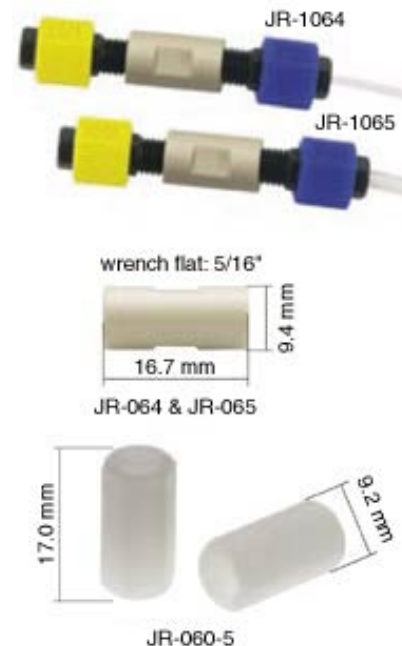

Jump to product numbers for:

[Unions](#)

[Reducing union](#)

[Spare parts](#)

TOP ▲

Low Pressure Unions

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8"


Description		Product No.
Polymide union, through bore, body only; 1/4"-28, for flanged tubing only, pkg of 5		JR-060-5
0.50 mm bore	PEEK union with fittings	JR-1064

	PEEK union, body only	JR-064
0.75 mm bore	PEEK union with fittings	JR-1068
	PEEK union, body only	JR-068
1.30 mm bore	PEEK union with fittings	JR-1065
	PEEK union, body only	JR-065

TOP ▲

PEEK Reducing Union

Description	Product No.
1/4" to 1/8", 1.5 mm bore, with fittings	JR-CU4LPK

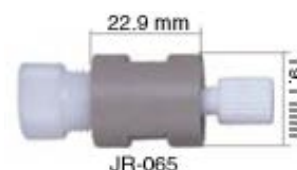

PEEK Union Spare Parts

Description		Qty/pkg	Product No.
Nuts	PPS, flangeless, 1/16", black	10	JR-55070-10
	PPS, flangeless, 1/8", black	10	JR-55071-10
	Polyacetal, flangeless, 1/8", white	10	JR-050WH-10
	Polyacetal, flangeless, 1/4", white	1	JR-CFL-4D
Ferrules	ETFE, 1/16"	10	JR-041-10
	ETFE, 1/8"	10	JR-051-10
	CTFE, 1/4"	1	JR-CFL-CB4KF-S

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

CUSTOMER SUPPORT

Tees, Crosses, and Manifolds

- Completely inert and biocompatible
- Different bore sizes available
- 1/4"-28 ports for 1/16" and 1/8" OD tubing

VICI tees, crosses, and manifolds for low pressure applications are machined from PEEK. The bodies have 1/4"-28 female flat bottom ports which accept any flanged fitting or collapsible or inverted flangeless fitting for 1/8" or 1/16" OD tubing.

Complete versions with 0.50, 0.75, and 0.80 mm bore are supplied with PPS flangeless nuts and ETFE ferrules for 1/16" OD tubing. Complete versions with 1.50 mm bore are supplied with PPS flangeless nuts and ETFE ferrules for 1/8" OD tubing. The fingertight sleeves are polypropylene.



Jump to product numbers for:

[Tees](#)

[Crosses](#)

[Manifolds](#)

[Spare parts](#)

Dimensions

Tees and Crosses

Body diameter: 25.0 mm; thickness: 12.0 mm

Manifolds

5 ports body diameter: 22.2 mm; thickness: 17.1 mm

9 ports body diameter: 22.2 mm; thickness: 28.6 mm

TOP ▲

Low Pressure Tees

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8".

Varies with tubing material and ID

Description		Product No.
0.50 mm bore	PEEK tee with 1/16" fittings	JR-6029
	PEEK tee, body only	JR-6034
0.80 mm bore	PEEK tee with 1/16" fittings	JR-6030
	PEEK tee, body only	JR-6032
1.50 mm bore	PEEK tee with 1/8" fittings	JR-6031
	PEEK tee, body only	JR-6033

TOP ▲

Low Pressure Crosses

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8".

Varies with tubing material and ID

Description		Product No.
0.50 mm bore	PEEK cross with 1/16" fittings	JR-6039
	PEEK cross, body only	JR-6044
0.80 mm bore	PEEK cross with 1/16" fittings	JR-6040
	PEEK cross, body only	JR-6042
1.50 mm bore	PEEK cross with 1/8" fittings	JR-6041
	PEEK cross, body only	JR-6043

TOP ▲

Low Pressure Manifolds

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8".

Varies with tubing material and ID

Description		Product No.
0.75 mm bore	5 port PEEK manifold with 1/16" fittings	JR-C5M1PK
	9 port PEEK manifold with 1/16" fittings	JR-C9M1PK

1.50 mm bore	5 port PEEK manifold with 1/8" fittings	JR-C5M2PK
	9 port PEEK manifold with 1/8" fittings	JR-C9M2PK

Spare Parts

Description		Qty/pkg	Product No.
Nuts	PPS, flangeless 1/4"-28, 1/16"	10	JR-55070-10
	PPS, flangeless 1/4"-28, 1/8"	10	JR-55071-10
Ferrules	ETFE, 1/16"	10	JR-041-10
	ETFE, 1/8"	10	JR-051-10

Also available: [Fingertight adapters in additional colors](#)

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

CUSTOMER SUPPORT

Color-Sleeve Fingertight Adapters

- Available in twelve colors
- Turns our PPS low pressure nut into a fingertight fitting


This sleeve is useful for color-coding different solvent lines, and turns the new VICI Jour low pressure nut into a fingertight fitting in an instant. The Color-Sleeve is available in 12 colors.

Color-Sleeve Fingertight Adapters

Material: Polypropylene

Color	Qty/ pkg	Product No.
Natural	10	JR-55080-10
White	10	JR-55081-10
Dark grey	10	JR-55082-10
Black	10	JR-55083-10
Lavender	10	JR-55084-10
Red	10	JR-55085-10
Yellow	10	JR-55086-10
Orange	10	JR-55087-10
Brown	10	JR-55088-10



Green	10	JR-55089-10	
Blue	10	JR-55090-10	
Purple	10	JR-55091-10	
2 each	24	JR-5508X-24	
2 blue 2 red 2 green 2 yellow 2 black 2 white	12	JR-5508X-12	

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops


Accessories/miscellaneous

CUSTOMER SUPPORT

Manual Flow Selection Valves

- Completely inert and biocompatible
- Several flow patterns available
- 8 mm bore, for use with 1/16" or 1/8" OD tubing
- Low dead volume
- For pressure up to 35 bar (500 psi)

VICI Jour manual valves are the ideal tools for liquid flow selection. The mountable valves are pre-fitted with the new VICI Jour PPS 1/4"-28 fittings and ETFE ferrules for 1/16" tubing, but the 1/4"-28 female flat bottom ports accept any flanged fitting or collapsible or inverted flangeless fitting for 1/8" or 1/16" OD tubing. Wetted parts are PTFE and CTFE; the body is polyacetal.

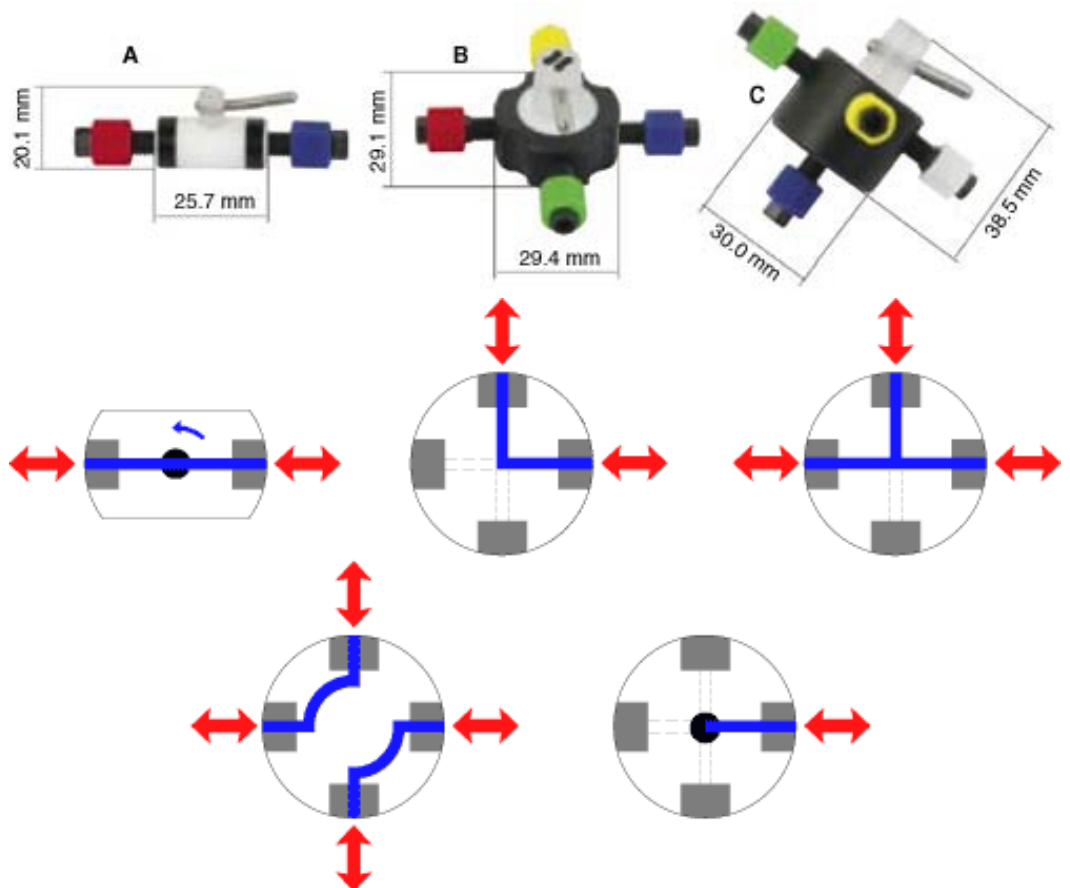

TOP ▲

Manual Flow Selection Valves

Pressure rating: < 35 bar (500 psi). Varies with tubing material and ID

No. of ports	Function	Photo	Product No.
2 ports	On/off	A	JR-MSV25
4 ports	L-flow	B	JR-MSV4L
	T-flow		JR-MSV4T
	2 pairs		JR-MSV4P
5 ports	1 inlet/4 outlets	C	JR-MSV5D

Spare Parts

Description		Qty/pkg	Product No.
Nuts	PPS, flangeless 1/4"-28, 1/16"	10	JR-55070-10
	PPS, flangeless 1/4"-28, 1/8"	10	JR-55071-10
Ferrules	ETFE, 1/16"	10	JR-041-10
	ETFE, 1/8"	10	JR-051-10

Also available: [Fingertight adapters in additional colors](#)

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

PPS Nuts for use with Collapsible Ferrules*


- Economically priced alternative to PEEK
- Excellent chemical resistance
- For 1/16" or 1/8" tubing, 1/4"-28 threads
- Slim design – fingertight with adapter

The new VICI Jour PPS 1/4"-28 fitting series offers excellent chemical and pressure resistance. The collapsible ferrule makes a very effective seal on a narrow area with no distortion of the tubing. No flanging tool is required.

*patent no. 6,575,501



1/8"

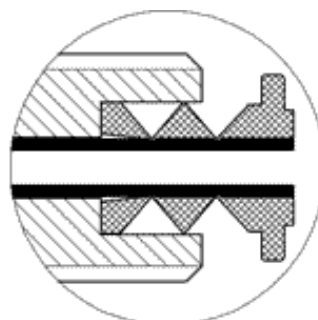

1/16"

[Jump to product numbers](#)

Chemical Resistance

Chemical Class	Resistance
Acids, organic	A
Acids, inorganic	A/C*
Aldehydes	A
Alcohols	A
Bases	A/B

* conc. Halogen acids

A = Suitable

B = Marginal; dependent on application

C = Not recommended

Esters	A
Halogenated org.	A/B
Hydrocarbons	A
Ketones	A

TOP ▲

PPS Nuts and Collapsible Ferrules

Pressure rating: 1/16" up to 105 bar (1500 psi); 1/8" up to 35 bar (500 psi). Varies with tubing material and ID

Material: PolyPhenylene Sulfide nut, CTFE ferrule

Description	Qty/ pkg	Product No.
1/16" nut	10	JR-55050-10
1/16" collapsible ferrule	10	JR-CFL-CB1KF
1/8" nut	10	JR-55051-10
1/8" collapsible ferrule	10	JR-CFL-CB2KF

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous


CUSTOMER SUPPORT

PPS Nuts for use with Flange/Washer Connection

- Economically priced alternative to PEEK
- Excellent chemical resistance
- For 1/16" or 1/8" tubing, 1/4"-28 threads
- Slim design – fingertight with adapter

The new VICI Jour PPS 1/4"-28 fitting series offers excellent chemical and pressure resistance. With this type of connection, the solvent contacts only the tubing.

For forming the flange on the tubing, we recommend our [Easy-Flange kit](#).

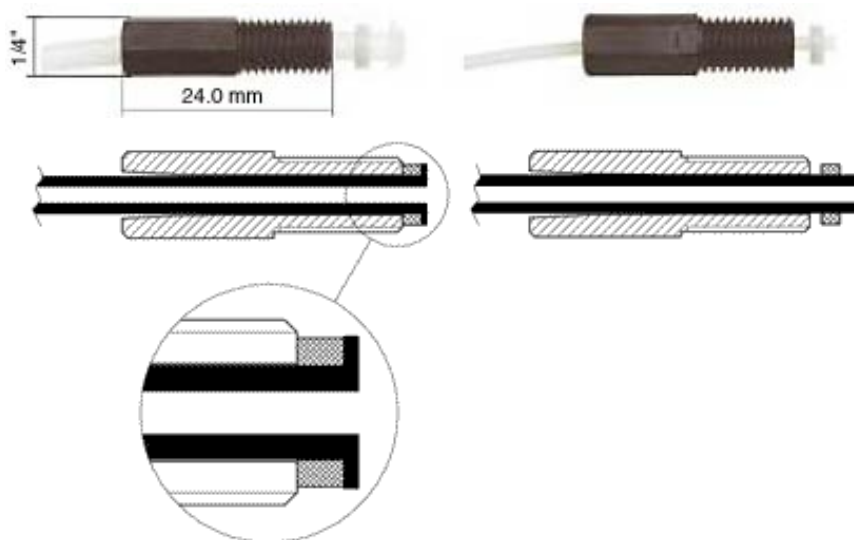

[Jump to product numbers](#)

Chemical Resistance

Chemical Class	Resistance
Acids, organic	A
Acids, inorganic	A/C*
Aldehydes	A
Alcohols	A
Bases	A/B
Esters	A

* conc. Halogen acids

A = Suitable

B = Marginal; dependent on application

C = Not recommended

Halogenated org.	A/B
Hydrocarbons	A
Ketones	A

TOP ▲

PPS Nuts and PP Washers for Flanged Connections

Pressure rating: up to 30 bar (435 psi). Varies with tubing material and ID

Material: PolyPhenylene Sulfide nut, Polypropylene washer

Description	Qty/ pkg	Product No.
1/16" nut	10	JR-55060-10
1/16" washer	10	JR-201926-10
1/8" nut	10	JR-55061-10
1/8" washer	10	JR-201928-10

TOP ▲

VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous


CUSTOMER SUPPORT

Polypropylene Color-Coded Flanged Nuts

- For flanged connections
- Compatible with all 1/4"-28 flat bottom ports
- Fiberglass reinforced polypropylene with polypropylene washers
- Available for 1/16" and 1/8" tubing

A number of chromatography systems still require or prefer flanged connections. When the tubing is flanged, the washer never comes in contact with the mobile phase. Six colors are available.

For forming the flange on the tubing, we recommend our [Easy-Flange kit](#).

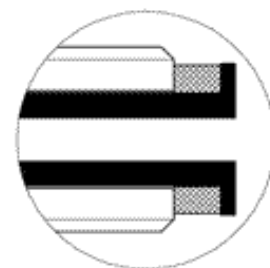

1/16" Polypropylene Color-Coded Flanged Nuts

Pressure rating: up to 30 bar (435 psi). Varies with tubing material and ID.

Description	Color	Qty/ pkg	Product No.
Polypropylene 1/4"-28 nut	Black	10	JR-201094BK-10
	White	10	JR-201094WH-10
	Blue	10	JR-202122-10
	Green	10	JR-202123-10
	Yellow	10	JR-202124-10
	Red	10	JR-202125-10
Polypropylene 1/4"-28 nut, 1" long	Black	10	JR-202118-10
	White	10	JR-202118-10

Washer for 1/16" flanged fitting	10	JR-201926-10
----------------------------------	----	--------------

TOP ▲

1/8" Polypropylene Color-Coded Flanged Nuts

Pressure rating: up to 30 bar (435 psi). Varies with tubing material and ID.

Description	Color	Qty/ pkg	Product No.
Polypropylene 1/4"-28 nut	Black	10	JR-201580BK-10
	White	10	JR-201580WH-10
	Blue	10	JR-202142-10
	Green	10	JR-202143-10
	Yellow	10	JR-202144-10
	Red	10	JR-202145-10
Polypropylene 1/4"-28 nut, 1" long	Black	10	JR-202119-10
	White	10	JR-202121WH-10
Washer for 1/8" flanged fitting		10	JR-201928-10

TOP ▲

VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

Polyacetal Color-Coded Flangeless Nuts

- Value-for-money alternative to PEEK
- Flange-free – compatible with all 1/4"-28 flat bottom ports
- Excellent chemical resistance
- Available for 1/16" and 1/8" tubing
- Fingertight – no tools required

Polyacetal flangeless nuts in combination with the inverted ferrule make leak-free seals without a flanging tool or a wrench. The inert ETFE ferrule is the only part besides the tubing which contacts the fluid – the nut itself does not contact the fluid stream. (See illustration) Six colors are available.



1/16" Polyacetal Color-Coded Flangeless Nuts

Pressure rating: up to 70 bar (1000 psi). Varies with tubing material and ID.

Description	Color	Qty/ pkg	Product No.
Polyacetal 1/4"-28 flangeless nut	Green	10	JR-040GR-10
	Blue	10	JR-040BL-10
	Yellow	10	JR-040YL-10
	Red	10	JR-040RE-10
	White	10	JR-040WH-10
	Black	10	JR-040BK-10
Polyacetal 1/4"-28 flangeless nut, 1" long	White	10	JR-040WL-10

1/16" ETFE ferrule	10	JR-041-10
--------------------	----	-----------

TOP ▲

1/8" Polyacetal Color-Coded Flangeless Nuts

Pressure rating: up to 35 bar (500 psi). Varies with tubing material and ID.

Description	Color	Qty/ pkg	Product No.
Polyacetal 1/4"-28 flangeless nut	Green	10	JR-050GR-10
	Blue	10	JR-050BL-10
	Yellow	10	JR-050YL-10
	Red	10	JR-050RE-10
	White	10	JR-050WH-10
	Black	10	JR-050BK-10
Polyacetal 1/4"-28 flangeless nut, 1" long	White	10	JR-050WL-10
1/8" ETFE ferrule		10	JR-051-10

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous


CUSTOMER SUPPORT

Flanged Tubing Connector

- For easy connection of flanged tubing
- Zero dead volume
- Available for 1/16" and 1/8" tubing

With its female 1/4"-28 thread, the VICI flanged tubing connector fits all of our flanged fittings. We recommend this PEEK connector instead of a union to make a zero dead volume connection.

For forming the flange on the tubing, we recommend our [Easy-Flange kit](#).

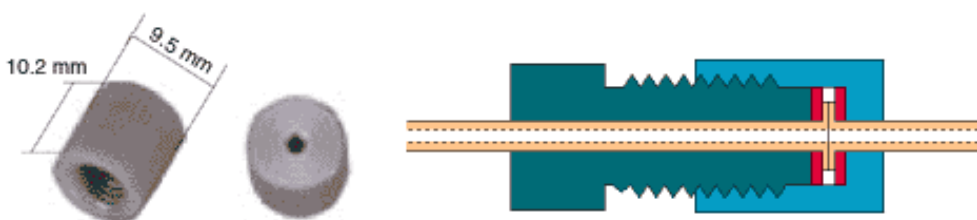

Flanged Tubing Connector

Pressure rating: up to 30 bar (435 psi). Varies with tubing material and ID.

Description	Product No.
For 1/16" tubing	JR-CEN1PK
For 1/8" tubing	JR-CEN2PK

VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

Low Pressure Bulkhead Unions

- Panel mounting
- Easy connection of low pressure solvent lines
- Low dead volume
- Completely inert and biocompatible



Low Pressure Bulkhead Unions

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8"

Description	Product No.
PEEK, 0.75 mm bore, with ETFE ferrules and PPS flangeless nuts for 1/16" tubing	JR- CBUMPK
PEEK, 1.50 mm bore, with ETFE ferrules and PPS flangeless nuts for 1/8" tubing	JR-CBULPK

Spare Parts

Description	Qty/pkg	Product No.
Nuts		
PPS, flangeless, 1/16", black	10	JR-55070-10
PPS, flangeless, 1/8", black	10	JR-55071-10
Ferrules		
ETFE, 1/16"	10	JR-041-10
ETFE, 1/8"	10	JR-051-10
Adapters		
PP fingertight sleeve, yellow	10	JR-55086-10
PP fingertight sleeve, blue	10	JR-55090-10



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

- Tubing for chromatography
- High pressure fittings
 - Nanovolume
 - Stainless steel
 - Polymeric
- Low pressure fittings
- Filters & mobile phase filters
- Mobile phase handling systems
- Frits
- PEEK columns & guard columns
- Sample loops
- Accessories/miscellaneous

CUSTOMER SUPPORT

JR-051-10

10

ETFE, 1/8"

Also available: Fingertight adapters in additional colors

PEEK Biocompatible Mixing Tees

- Completely inert and biocompatible
- Mountable
- Two different bore sizes available
- 1/4"-28 ports for 1/16" and 1/8" OD tubing



A unique mixing is the result of the turbulence generated by the introduction angle of the solvents. The streams mix more thoroughly than with a conventional tee.

The 0.75 mm bore tee is supplied with PPS flangeless nuts and ETFE ferrules for 1/16" OD tubing. The 1.50 mm bore version come with PPS flangeless nuts and ETFE ferrules for 1/8" OD tubing. The fingertight sleeves are polypropylene.

Dimensions

Body length: 44.5 mm; thickness: 22.2 mm

PEEK Biocompatible Mixing Tees

Pressure rating: up to 105 bar (1500 psi) for 1/16", 35 bar (500 psi) for 1/8".

Varies with tubing material and ID

Description	Product No.
PEEK mixing tee, 0.75 mm bore, with 1/16" fittings	JR-CM1XPK
PEEK mixing tee, 1.50 mm bore, with 1/8" fittings	JR-CM2XPK

Spare Parts

Description		Qty/pkg	Product No.
Nuts	PPS, flangeless 1/4"-28, 1/16"	10	JR-55070-10
	PPS, flangeless 1/4"-28, 1/8"	10	JR-55071-10
Ferrules	ETFE, 1/16"	10	JR-041-10



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

Low Pressure Polypropylene Adapters

- Adapt 10-32, 1/4"-28, barb, Luer, etc
- Completely inert and biocompatible

Description		Qty/ pkg	Product No.
Male Luer to:	1/16" barbed	10	JR-070211-10
	3/32" barbed	10	JR-070212-10
	1/8" barbed	10	JR-070213-10
	10-32 male	10	JR-070224-10
	14"-28 female	10	JR-070225-10
Female Luer to:	1/16" barbed	10	JR-070111-10
	3/32" barbed	10	JR-070112-10
	1/8" barbed	10	JR-070113-10
	Female Luer	10	JR-070101-10
	10-32 male	10	JR-070124-10
	1/4"-28 male	10	JR-070123-10
1/4"-28 male to:	1/16" barbed	10	JR-071123-10
	1/8" barbed	10	JR-071323-10



JR-071124-10	JR-071324-10	JR-072450-10	JR-072350-10
10	10	10	10
1/16" barbed	1/8" barbed	10-32	1/4"-28
10-32 male to:	Plugs		



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

Tubing for chromatography

High pressure fittings

- Nanovolume
- Stainless steel
- Polymeric

Low pressure fittings

Filters & mobile phase filters

Mobile phase handling systems

Frits

PEEK columns & guard columns

Sample loops

Accessories/miscellaneous

CUSTOMER SUPPORT

Peristaltic Tubing Adapter Set

- Connects soft tubing to PTFE and PTFE-like tubing
- Works with peristaltic tubing up to 4.50 mm OD
- Biocompatible
- Low dead volume
- For 1/16" and 1/8" OD PTFE and PTFE-like tubing

This "one-piece-fits-all" set consists of everything needed to connect 1/16" and 1/8" OD PTFE and PTFE-like tubing to soft tubing like silicone up to 4.50 mm OD. The adapter comes with three different bushings (for the different soft tubing OD's) and two PPS fittings with inverted ETFE ferrules (one each for 1/16" and 1/8" OD tubing). The body and bushings are PEEK.

Dimensions

Length (without fitting): 28 mm; diameter: 12 mm; bore: 0.50 mm



Peristaltic Tubing Adapter Set

Description	Product No.
PEEK adapter for peristaltic to PTFE tubing connection	JR-PTAS

TOP ▲



VICI JOUR - Low Pressure Fittings

VICI JOUR PRODUCTS

[Tubing for chromatography](#)

[High pressure fittings](#)

- [Nanovolume](#)
- [Stainless steel](#)
- [Polymeric](#)

[Low pressure fittings](#)

[Filters & mobile phase filters](#)

[Mobile phase handling systems](#)

[Frits](#)

[PEEK columns & guard columns](#)

[Sample loops](#)

[Accessories/miscellaneous](#)

Flanged Fittings Kit

- A complete kit of flanged fittings and accessories
- All parts biocompatible
- Excellent chemical resistance

The flanged fitting kit contains all the materials that are required to make flanged connections for 1/16" and 1/8" OD tubing: flanged fittings, washers, unions, tubing in four different sizes, our Clean-Cut tubing cutter, and a pair of tweezers.



CUSTOMER SUPPORT

Flanged Fittings Kit

Description	Product No.
Flanged fittings kit	JR-201529

[Complete list of contents \(pdf\)](#)

[Home](#) | [Products](#) | [Search](#) | [Contact VICI Jour](#) | [About VICI Jour](#)

VICI® is a registered trademark of Valco Instruments Co. Inc. and VICI AG
© 2005 VICI

Index of contents of the Flanged Fitting Kit (Part-No. JR-201529)

Part-No.	Description	Qty
JR-201094Black	Black flanged nut 1/16"	10
JR-201094White	White flanged nut 1/16"	10
JR-201580Black	Black flanged nut 1/8"	10
JR-201580White	White flanged nut 1/8"	10
JR-201928Bushing	1/8" Polypropylene bushing	20
JR-201926Bushing	1/16" Polypropylene bushing	20
JR-060	Union	5
JR-797	Clean-cut tubing cutter	1
JR-4011	1/16" OD x 0.25 mm ID PTFE tubing	3 m
JR-4183	1/16" OD x 0.50 mm ID PTFE tubing	3 m
JR-4036	1/16" OD x 0.75 mm ID PTFE tubing	3 m
JR-6800	1/8" OD x 1.59 mm ID PTFE tubing	3 m
JR-806	Tweezer tool	1