

## Selectors – Low Pressure

### Stream selector, 1/16" Valco ZDV fittings, 0.75 mm ports (.030")

Model C35Z

**Low pressure**

Includes Valco ZDV  
PEEK nuts and ferrules.


**Stream selector**

**10-32 ZDV**

Available only with microelectric actuator:  
24 VDC, with 110/230 VAC to 24 VDC  
power supply.

**1/16" 0.75 mm**

(For clarity, only eight of  
the twenty-four streams  
are illustrated.)

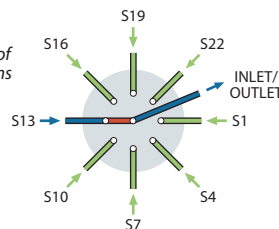


#### SPECS

**100 psi liq**  
**50°C max**  
PPS stator  
Valcon E2 rotor

#### OPTIONS

- Optional bore:  
0.5 mm (.020")  
1.0 mm (.040")
- 2", 3", 4", and 6"  
standoffs
- Consult the factory  
for optional materials.

	20 Position		24 Position		26 Position	
	Prod No	Price	Prod No	Price	Prod No	Price
With microelectric actuator	C35Z-31820EMT		C35Z-31824EMT		C35Z-31826EMT	
Replacement valve	C35Z-31820D		C35Z-31824D		C35Z-31826D	
Replacement rotor	C35Z-31R20		C35Z-31R20		C35Z-31R20	

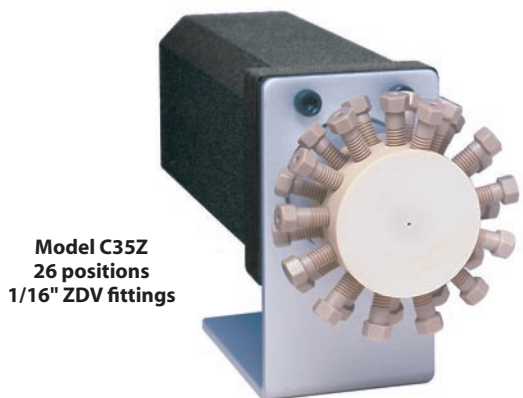

**Model C35Z**  
**26 positions**  
**1/16" ZDV fittings**


#### OPTIONAL FLOWPATHS

**Model C35Z** valves select and isolate one of 20-26 streams, with the remainder dead-ended.


**Model C35ZF**, the flow-through version, is similar to the C35Z but its non-selected streams continue flowing through individual outlets. 10, 12, and 13 positions are available.

**Model C35ZT**, the trapping version, is similar to the C35ZF but has a second selected port. Non-selected streams continue flowing. 10, 12, and 13 positions are available.

Call for pricing and information.

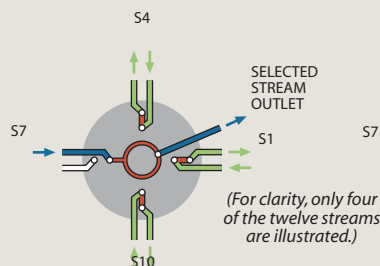

**Model C35ZF**  
**schematic**

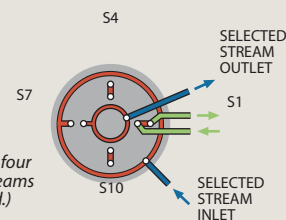

**Model C35ZT**  
**schematic**