

New from VICI

NEW

360 Micron Fittings pages 57-58

Our new high pressure fittings permit direct connection of 360 micron OD fused silica, PEEK, stainless, or electro-formed nickel tubing without the use of liners. The ferrule snaps into the nut so that the fitting is "one-piece", but the ferrule remains free to rotate as the nut is tightened so that the tube doesn't twist. Because of the compact size and fine 2-56 threads, a leak-free connection that seals at pressures in excess of 20,000 psi can be easily formed with the available manual tool.




360 micron nanovolume fittings
page 57-58

NEW

Chem inert UHPLC Injectors and Selectors Model C70 series

- Nanovolume®**pages 152-153
- Microbore** 156-157
- Selectors** 170-171

The ideal injector for current analytical techniques requiring system pressures up to 20,000 psi, C70 series injectors come in 4, 6, 8, and 10 port versions, with bores as small as 100 µm (.004"). Selectors offer 4 to 10 positions. In addition to injectors incorporating our new 360 µm fittings, models are available for use with 1/16" or 1/32" tubing.

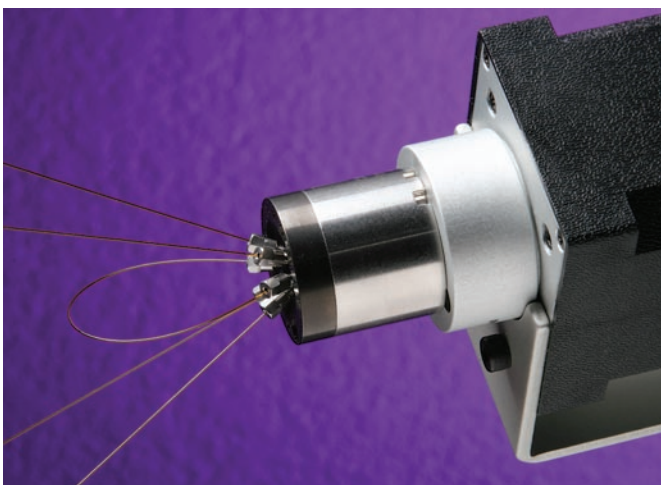

UHPLC injector with 360 micron nanovolume fittings
Model C72MU, page 152