



Fused Silica Adapters

Fused silica make-up adapters

The fused silica make-up adapter connects a fused silica capillary column to a valve or detector while adding a make-up gas. In the reverse mode it works like a splitter, without the uneven or erratic split seen with basic tees. Two lengths are available. Order 1/32" fused silica adapter ferrules separately (see box on facing page).

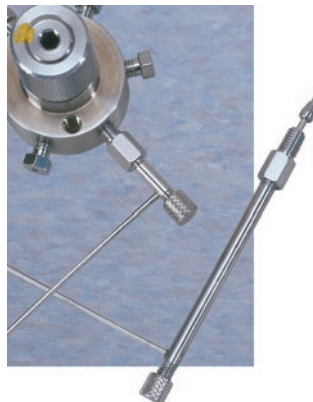

Description	Length	Bore	Prod No	Price
1/16" to 1/32"	1.5"	0.5 mm	FSMUAS1.5M	
	1.5"	0.75 mm	FSMUAS1.5	
	1.5"	1.0 mm	FSMUAS1.5L	
	3.5"	0.75 mm	FSMUA1.5	

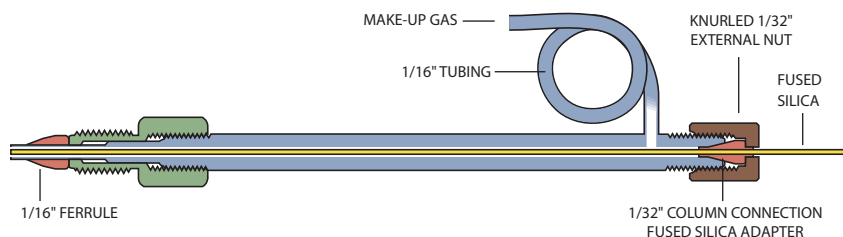

Fused silica make-up adapter
(FSMUA1.5)

100 μm	= .004"
150 μm	= .006"
0.25 mm	= .010"
0.50 mm	= .020"
0.75 mm	= .030"
1.0 mm	= .040"
1.5 mm	= .060"
2.0 mm	= .080"
4.6 mm	= .180"
6.0 mm	= .236"
6.4 mm	= .253"
7.0 mm	= .275"
10.0 mm	= .400"
27.0 mm	= 1.08"
1/32"	= 0.8 mm
1/16"	= 1.6 mm
1/8"	= 3.2 mm
1/4"	= 6.4 mm
3/8"	
	= 9.5 mm
1/2"	= 12.7 mm