

Restek... Leading the Pack in Protection




How many times has it been said about press-fit connectors, "Some people swear by them ... others swear at them"? For many analysts, the art of attaching a guard column to an analytical column is a mystery. Restek's chemists have discovered the solution to this mystery—the most reliable connection is no connection at all! No guard column system is more permanent than one continuous length of tubing containing both the guard column and the analytical column.

Restek now offers a wide variety of *Integra-Guard*™ capillary columns. A guaranteed leak-free connection! The column is tied using high temperature string with the guard

column set apart from the analytical column. The transition area between the columns is the point at which the guard column ends and the analytical column begins. The entire setup is suspended in our unique "crush-free" cage that prevents the column from coming in contact with anything that could damage it.

The *Integra-Guard*™ column is so economical that we challenge you to beat our price with that of a conventional connection, even if you assemble it yourself. If you are currently using a guard column or considering one for the future, call Restek today and ask about our *Integra-Guard*™ columns.

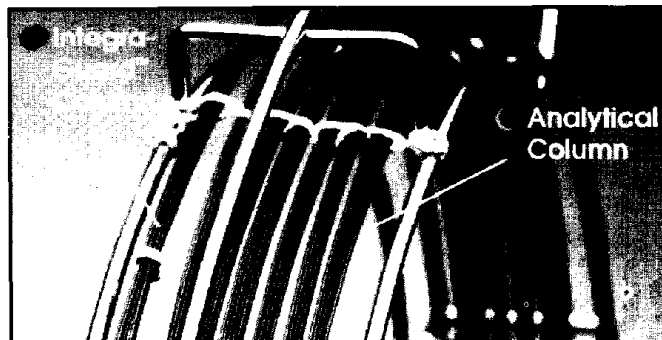

Ordering is easy. Just add the appropriate suffix number to the analytical column's catalog number. For example, a 30m, 0.25mmID, 0.25um Rtx™-5 with a 5meter *Integra-Guard* Column is cat.# 10223-124.

mmID	Length	Suffix #
0.25	5m	-124
	10m	-127
0.32	5m	-125
	10m	-128
0.53	5m	-126
	10m	-129

Protecting your analytical column has never been this easy!

INTEGRA-GUARD

INTEGRATED GUARD COLUMNS

